

Databrowse: An extensible data management platform

Tyler J. Lesthaeghe
tylerl@iastate.edu

Nathan D. Scheirer
scheirer@iastate.edu

Department of Aerospace Engineering and
Center for Nondestructive Evaluation
Iowa State University of Science and Technology
Ames, Iowa 50011

June 19, 2018
Version 0.8

Databrowse: An Exensible Data Management Platform

Copyright © 2012-2016 Iowa State University Research Foundation, Inc.

All rights reserved.

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

1. Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.
3. Neither the name of the copyright holder nor the names of its contributors may be used to endorse or promote products derived from this software without specific prior written permission.

THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AND CONTRIBUTORS “AS IS” AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT HOLDER OR CONTRIBUTORS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

This material is based on work supported by the Air Force Research Laboratory under Contract #FA8650-10-D-5210, Task Order #023, and performed at Iowa State University.

DISTRIBUTION A. Approved for public release: distribution unlimited; 19 Aug 2016; 88ABW-2016-4051.

This material is based on work supported by NASA under Contract NNX16CL31C and performed by Iowa State University as a subcontractor to TRI Austin.

Approved for public release by TRI Austin: distribution unlimited; 01 June 2018; by Carl W. Magnuson (NDE Division Director).

Contents

1	Introduction	5
1.1	Motivation	5
1.2	Implementation	5
1.2.1	XML	5
1.2.2	XSLT	6
1.3	Databrowse Plugins	6
2	Installing Databrowse	9
2.1	Databrowse Platforms	9
2.2	Databrowse Server	9
2.2.1	Software Requirements and Dependencies	9
2.2.2	Installation on Unix-based Platforms	10
2.3	CEFDatabrowse	12
2.3.1	Software Requirements and Dependencies	12
2.3.2	Installation	13
3	Getting Started	14
3.1	Features of Interest	14
3.2	Creating New Directories	15
3.3	Uploading Files	16
4	Databrowse Configuration Options	17
4.1	databrowse.wsgi.conf	17
4.2	iconmap.conf	18
4.3	hiddenfiles.conf	18
5	Databrowse Library	19
6	Databrowse Plugin Format	20
6.1	Plugin File Structure	20
6.2	File Contents	20
6.2.1	__init.py__	20
6.2.2	db_plugin_name.py	20
6.2.3	dbs_stylesheet_name.xml	23
6.2.4	handlers.py	24
7	Included Databrowse Plugins	25
7.1	Checklist Editor	25
7.2	Checklist Viewer	25
7.3	Datacollect v1 Viewer	26
7.4	Datacollect v2 Viewer	26
7.5	Dataguzzler Data File	26
7.6	Dataguzzler Settings File	27
7.7	Data Table	27
7.8	Default	27
7.9	Directory	28
7.10	File Operations	28
7.11	Generic Binary File	28
7.12	Generic HDF5 File	29
7.13	Generic WSGI Application	29
7.14	Generic XML File	29
7.15	Image Viewer	30
7.16	Mercurial Repository	30

7.17	Movie Viewer	30
7.18	Multimedia Directory	31
7.19	Office Viewer	31
7.20	PDF Viewer	32
7.21	Plain Text File	32
7.22	SolidWorks Viewer	32
7.23	Specimen Management Plugin	33
7.24	Specimen Directory Plugin	33
7.25	Specimen Group Management Plugin	33
7.26	SVG Viewer	34
7.27	Transducer Management Plugin	34
7.28	Transducer Directory Plugin	34
7.29	Trigger Log Plugin	34
7.30	Trigger Log Directory Plugin	34
7.31	Web Page Viewer	35
8	License and Third-Party Components	36
8.1	BSD 3-Clause License	37
8.2	GNU General Public License v3	38
8.3	GNU Lesser General Public License v3	47
8.4	MIT License	50
8.5	Zope Public License v2.1	51
8.6	Apache License v2.0	52

1 Introduction

Databrowse is an extensible web-based platform for data viewing, manipulation, and management. At the most basic level, Databrowse is a file browser; however, simple plugins enable Databrowse to represent data of a variety of formats in a consistent way, enabling rapid viewing and transformation of those views to narrow in on features of interest. The plugin architecture enables Databrowse to be adapted to support any data format in which knowledge of the data format is available. Data from multiple sources or formats can be pulled together into combined representations and then further transformed as desired. Furthermore, all of these transformations are performed in real time.

1.1 Motivation

Databrowse was originally developed to aid in viewing and analyzing data collected in the field of nondestructive evaluation (NDE). NDE is a broad, highly interdisciplinary field related to the development of measurement techniques that find and characterize material flaws and condition. Some well known NDE techniques include visual/liquid penetrant inspection, magnetic particle inspection, ultrasonics, radiography, and eddy current testing.

NDE techniques are capable of generating considerable amounts of data in very short periods of time. However, many industrial NDE inspections today produce a simple pass/fail response as a result of a testing process. Recorded raw data is often viewed as useless, or potentially even a liability, unless we have ways to extract useful information. Ultimately, the problem of finding a needle in a haystack is particularly challenging, especially if you do not know what you are looking for.

Even in research scenarios where we desire to collect large quantities of data to examine specific items, handling large quantities of data can still be a challenge. As a part of a recent modeling effort related to development of a forward model for vibrothermography, a nondestructive testing technique that utilizes vibration-induced heating to locate cracks in materials, the need became apparent for better tools and data management practices. It was known at the beginning of the work that a considerable amount of data was to be generated. In the end, a final data table containing almost 25,000 entries was generated, along with over 0.5 TB of raw data.

The authors sought to develop a tool that would help in the short term for dealing with the latter problem, while serving as a spring board to further work toward dealing with the former problem. Databrowse was the resulting tool.

1.2 Implementation

1.2.1 XML


Databrowse represents data as XML. XML (eXtensible Markup Language) is a standard that allows text data to be hierarchically structured utilizing arbitrarily defined tags. Special engines can then be used to parse and manipulate such structured data. Figure 1 shows an example of a simple data set being structured and stored in an XML format. The first line of the sample XML file shows an experiment tag, and the last line shows that tag being closed. Everything contained between the opening and closing tags can be described as children. In this context, our measurements are children of the experiment. Thus, a hierarchical data structure can be developed. Subsequently, we have represented all of the individual parameters for each measurement as children of that measurement.

We can take it a step further and indicate that one of these parameters could have multiple values, as seen with the voltage data. This is one major advantage of this type of data structure, as our spreadsheet style data table does not necessarily make representing this type of data structure easy. We have made it work here with a comma list of values in our data table; however, this does not work so easily with more complicated data.

Furthermore, XML tends to work very nicely as a way of representing most data, since many frequently used data formats internally represent data in a natural hierarchical structure, often as a convenient way of dealing with the issue just described. Therefore, conversion to XML, if even required, is generally very straightforward. Once represented as XML, Databrowse is able to leverage the power and speed of the open

#	Sample	Voltage	Result
1	A-001	1V, 2V, 3V	10
2	A-002	1.5V	15
3	B-004	1.5V, 3V	12

Simple Experimental Data Table

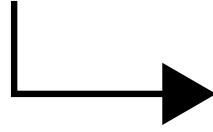

Converted and Stored in a
Simple XML Representation

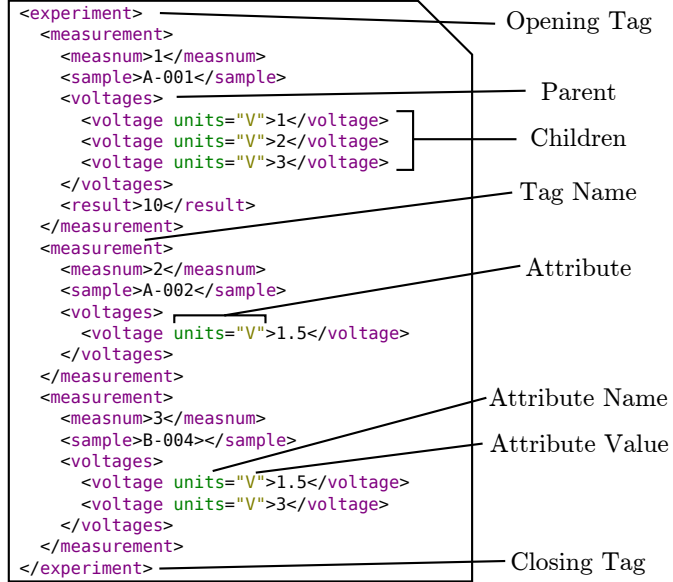


Figure 1: Sample XML Data

source XML engine libxml2. The result is being able to parse and transform considerable amounts of data in seconds.

It is important to note that raw binary waveforms are not intended to be represented as XML; however, pieces of them might be. More typically though, a plugin would utilize an interface provided by Databrowse that enables the creation of images of such data. Such an image can be generated in real time and served to the web browser by Databrowse. This would not be limited to images. Any format that could be displayed in a web browser could be used, such as videos, animations, or any format that can use a web browser plugin to display.

1.2.2 XSLT

Databrowse utilizes XSLT to transform data. XSLT (eXtensible Stylesheet Language Template) provides an interface by which the user can define a set of transformations to be applied to XML data. In other words, XML data and XSLT templates designed to act on that data are provided to the XSLT engine and the engine will output a new set of XML data based on the transformation provided. This behavior can be seen in Figure 2, where our sample data from Figure 1 has been transformed using an XSLT template.

Databrowse, being a web based platform, wants to build web pages utilizing the data being provided. HTML, the language with which web pages are built, is a type of XML. As a result, our XSLT transform is able to take our data file and build a web page dynamically in real time. Databrowse then handles the process of serving that web page to a user in their web browser.

1.3 Databrowse Plugins

Databrowse plugins are responsible for providing the following: 1) registering the file types with which they should be able to operate on, 2) providing an XML representation of the file, and 3) providing an XSLT transform that converts the XML representation to the desired HTML view that is displayed in the web browser. Plugins are also able to provide additional features that can be triggered or accessed from the web view. Such features might include the automatic generation of animations, data conversions, or running of processing scripts.

Since every file can be represented as XML, Databrowse provides an interface for recursively obtaining an XML representation of entire directories and sub-directories. Thus, a single representation of multiple files can be obtained. In addition to providing a set of plugins for some common file types, Databrowse

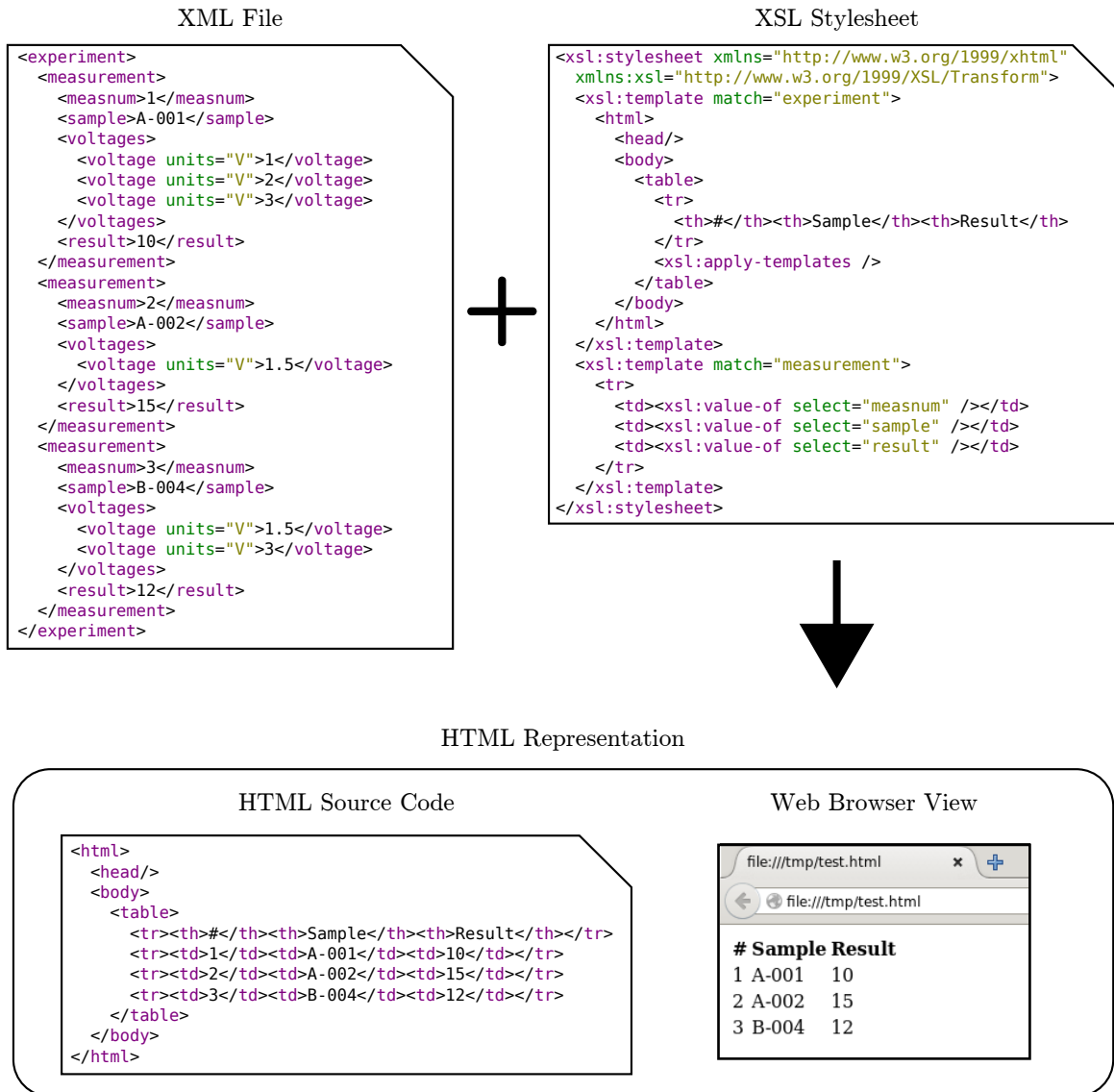


Figure 2: Sample XML Data being Transformed with XSLT

includes some plugins designed to provide a simple interface for building such combined representations. XSLT transformation stylesheets can also be provided on a per-directory basis if additional control is needed for specific use cases.

2 Installing Databrowse

2.1 Databrowse Platforms

Databrowse comes in two flavors, both of which contain the same functionality and provide nearly identical feature sets. Databrowse Server utilizes a WSGI compliant web server which is hosted the end users data storage system. The server version of Databrowse allows the end user to access their file system remotely from any other device connected to the internet following a onetime installation process. Alternatively, CEF-Databrowse has been implemented which utilizes CEFpython (Chromium Embedded Framework - Python) to provide the entirety of the Databrowse library on a single client. This self contained client can be installed directly through PIP and will seclude the Databrowse session to the host machine. This version of Databrowse can be especially effective in environments where data cannot be accessed via the internet or be even exposed to a network.

2.2 Databrowse Server

2.2.1 Software Requirements and Dependencies

Databrowse requires the use of a WSGI compliant web server. Databrowse has been tested extensively with mod_wsgi on Apache 2.2.15 on Red Hat Enterprise Linux 6.6 and Apache 2.4.10 on Ubuntu 15.04. The usage of mod_wsgi, Apache, and a Unix-based operating system are strongly recommended. Using Apache on Windows has been preliminarily tested; however, its use is strongly discouraged at this time due to file path issues that will be resolved in a later version of Databrowse.

Modern versions of Mozilla Firefox and Google Chrome are recommended for accessing Databrowse. Usage of Microsoft Internet Explorer is strongly discouraged. Microsoft Edge has not yet been tested.

Databrowse is dependent upon Python 2 (2.6 or later, Python 3 is not supported at this time). Databrowse also requires the following Python modules (available from PIP or the package management systems on Red Hat and Ubuntu):

- python-lxml version 3.2 or greater
- python-magic version 0.2 or greater (python-magic-bin if platform is Windows)
- python-numpy version 1.8 or greater
- python-pillow version 2.3 or greater

Several Databrowse plugins packaged with the distribution also require the use of the following Python modules:

- python-qrcode version 4.0 or greater
- Dataguzzler Python Bindings (<http://thermal.cnde.iastate.edu/dataguzzler>)
- Dataguzzler Units Support (<http://thermal.cnde.iastate.edu/dataguzzler>)

Databrowse also requires the use of mod_rewrite for URL rewriting. Databrowse can be ran without this support enabled; however, its use has not been tested extensively.

2.2.1.1 Experimental Windows Use A set of binary packages providing the necessary prerequisites have been compiled and are available upon request for usage with Apache on Windows. These packages require the usage of Apache compiled against the Microsoft Visual C++ 9 Runtime. This is a limitation imposed by Python 2. WampServer 2.2d is one such binary distribution of Apache that meets these requirements. Please contact the author for additional information.

2.2.2 Installation on Unix-based Platforms

The necessary minimum prerequisites can be installed using the following commands on Debian-based platforms:

```
sudo apt-get install python2.7 python-lxml python-magic (python-magic-bin if platform \
is Windows) python-numpy python-pip libjpeg libjpeg-dev libfreetype6 \
libfreetype6-dev zlib1g-dev apache2 libapache2-mod-wsgi libapache2-mod-rewrite \
python-dev python-setuptools libtiff5-dev libjpeg8-dev zlib1g-dev liblcms2-dev \
libwebp-dev tcl8.5-dev tk8.5-dev python-tk \
sudo pip install pillow
```

And for Fedora-based platforms:

```
sudo yum install python python-lxml python-magic python2-numpy python-pip libjpeg-turbo \
libjpeg-turbo-devel freetype freetype-devel zlib zlib-devel httpd mod_wsgi \
python-devel redhat-rpm-config libtiff-devel libjpeg-devel zlib-devel freetype-devel \
lcms2-devel libwebp-devel tcl-devel tk-devel
sudo pip install pillow
```

The package `redhat-rpm-config` is only required on Fedora 23 or later. The command `dnf` may also be exchanged for `yum` on newer Fedora-based platforms as well. Obtain the Databrowse source from the Databrowse website or from the GitHub repository (coming soon). From within the root directory, run the following command:

```
sudo python setup.py install
```

This will install the Databrowse library components into the the Python site-packages directory. Ensure that the necessary Apache modules are running and available by running the following command for Debian platforms:

```
sudo a2enmod wsgi
sudo a2enmod rewrite
```

And for Fedora platforms, both modules should be enabled by default. Verify the presence of the following lines in the file `/etc/httpd/conf/httpd.conf` or in any files contained in `/etc/httpd/conf.d/` or `/etc/httpd/conf.modules.d`:

```
LoadModule rewrite_module modules/mod_rewrite.so
LoadModule wsgi_module modules/mod_wsgi.so
```

From the Databrowse source folder, copy the contents of the `databrowse_wsgi` folder to an appropriate location. For example, to install the server components in `/var/www/databrowse`, use the following command from within the Databrowse source folder:

```
sudo mkdir /var/www/databrowse
sudo cp -a databrowse_wsgi/* /var/www/databrowse/
sudo chmod -R 755 /var/www/databrowse
```

Apache must now be configured. If working on a Debian-based system with a newer version of Apache, create the file `/etc/apache2/sites-available/databrowse.conf` with the following example contents and adjust as necessary for your system. On Fedora-based systems, the file name should be `/etc/httpd/conf.d/databrowse.conf` with the contents below adjusted as necessary:

```
WSGIScriptAlias /databrowse /var/www/databrowse/databrowse.wsgi
Alias /dbres /var/www/databrowse/resources
```

```

<Location "/databrowse">
    # Implementation of proper user controls is strongly encouraged!

    # Require SSL is also Strongly Encouraged But Must Be Appropriately Configured
    # Also be sure to modify databrowse_wsgi.conf to reflect 'https'
    # instead of 'http' on all URLs

    # SSLRequireSSL

    # Suggest the following items to restrict to localhost
    # Only disable if you know what you're doing and have Apache properly configured
    # AND have authentication enabled on Databrowse
    Order deny,allow
    Deny from all
    Allow from 127.0.0.1

    Options FollowSymLinks

    # Rewrite Rules - No Modification Should Be Needed Unless You Change Location Above
    RewriteEngine on
    RedirectMatch ^databrowse$ /databrowse/
    RewriteCond %{QUERY_STRING} path=
    RewriteRule ^(.*)/databrowse\.wsgi(.*)$ $1/databrowse.wsgi [QSA]
    RedirectMatch ^databrowse$ /databrowse/
    RewriteCond %{QUERY_STRING} !path=
    RewriteRule ^(.*)/databrowse\.wsgi(.*)$ $1/databrowse.wsgi?path=/$2 [QSA,L]
</Location>

```

WARNING! Improper configuration of a web server can leave your computer and data at risk of exposure. Databrowse does not provide any built-in authentication mechanism nor should it be relied on to prevent access to files outside of the configured data root. Accordingly, usage of Databrowse on computer systems in which web server access is available from the Internet is strongly discouraged without the usage of SSL and an appropriate Apache authorization module. It is strongly encouraged that anyone deploying Databrowse for use over the Internet be comfortable with securely configuring and using Apache. The authors provide no warranty and are not responsible for loss or damages. Please see the Databrowse license for additional information.

The above configuration assumes that the Databrowse WSGI components are installed in `/var/www/databrowse` and will serve Databrowse from `http://localhost/databrowse` and will serve Databrowse static resources from `http://localhost/dbres` on a default installation of Apache on Debian-based systems. The above configuration will utilize server default access permissions, so, you may wish to add additional configuration or ensure that your web server is not accessible from the Internet. Refer to the warning above.

The following command is required for Debian-based systems only and can be used to enable Databrowse:

```
sudo a2ensite databrowse.conf
```

Prior to restarting Apache, the files `databrowse_wsgi.conf` and `databrowse_style.xml` (or symbolic links to these files) must exist in the same directory as `databrowse.wsgi`. It should be sufficient to simply rename the file `databrowse_style.sample.xml` to `databrowse_style.xml` until you are ready to customize the appearance of Databrowse. The `databrowse_wsgi.conf` file will require some additional configuration. A sample file `databrowse_wsgi.sample.conf` is provided.

The most critical line that must be changed in `databrowse_wsgi.conf` from the sample provided file is `self.dataroot`. This should be set to the absolute path of the directory containing the data you wish Databrowse to serve. Databrowse provides very simple checks to keep the user inside directories below this path; however, this should not be used as a method of security! Symbolic links pointing out of this directory will be followed and there may be other methods available for a user to escape out of the data root path.

Additionally, you may wish to set `self.requireuser=False` if you do not presently have authentication configured on Apache. This may be okay if your Apache instance is only accessible from the local machine. However, if it is accessible to the world, it is highly encouraged to correctly configure authentication. See the previous warning about security.

If you have SSL properly configured on your web server and you uncommented the `SSLRequireSSL` line in `databrowse.conf`, you should update all URLs in `databrowse.wsgi.conf` to show `https` at the beginning instead of `http`.

It is also worth noting that files contained within data root must be at least readable by the web server process owner. To take advantage of the full capabilities of Databrowse, you will also want write access available to the web server process owner as well. Again, use caution when using Databrowse, especially if you are not familiar with securely configuring Apache. Additionally, if using SELinux or similar systems, you will need to make additional configuration changes to ensure the web server process has read/write access to all files in the data repository. On newer Fedora systems, for example, if you set your `dataroot` to `/home/databrowse`, then you could run the following command to allow access to Apache:

```
sudo chcon -R -t httpd_sys_rw_content_t /home/databrowse
```

Within `databrowse.wsgi.conf`, you will also want to ensure that `self.siteurl` and `self.resurl` are correct if you have changed any Apache settings from the defaults. You should also ensure that the directory listed in `self.checklistpath` has been created within the data root directory. This setting is related to the Checklist plugin, but, Databrowse presently will not run if this directory does not exist. This will be corrected in future versions of Databrowse.

Please see Chapter 4 for more information on these files.

Restart the web server and Databrowse should now be ready for use. On Debian-based systems:

```
sudo service apache2 restart
```

On Fedora-based systems:

```
sudo service httpd restart
```

You should now be able to access Databrowse by visiting `http://localhost/databrowse` in your web browser.

2.3 CEFDatabrowse

CEFDatabrowse is a standalone client that operates nearly identically to the server version of Databrowse, however this version does not require a connection to a network and allows local usage of Databrowse's feature set.

2.3.1 Software Requirements and Dependencies

CEFDatabrowse utilizes the same dependencies as the server version with a few additions. The Databrowse core architecture has been extensively tested on Linux, however the CEFDatabrowse client is being used on both Linux and Windows platforms at this time. Your experience may be different if you are using an untested platform. CEFDatabrowse requires Python 2.7 and will not function with any other version of Python at this time. The following Python modules are also required, these should be automatically installed when the PIP package or `setup.py` file is ran but if a package does not install correctly these are what are required:

- python-lxml version 3.2 or greater
- python-magic version 0.2 or greater (python-magic-bin if platform is Windows)
- python-numpy version 1.8 or greater
- python-pillow version 2.3 or greater
- cefpython3 version 57.0 or greater

2.3.2 Installation

On either Linux or Windows the installation process should be almost identical with the exception that the python-magic library should be used on Linux and the python-magic-bin library should be used on Windows. To install CEFDatabrowse verify that Python 2.7 is installed on your system. If not, run:

Linux:

```
sudo apt-get install python2.7
```

Windows:

Install Anaconda package manager or a Python 2.7 binary

CEFDatabrowse can be installed via PIP or from source. The PIP binaries are verified to work in both Linux and Windows environments, therefore to install via PIP execute the following command in your command prompt:

```
pip install databrowse
```

This will install CEFDatabrowse along with the components that would be required to run the server version to your sites-packages directory within your Python installation directory. Following a successful installation CEFDatabrowse can be utilized in the following manner:

Usage: databrowse [-h] [-s path] [-e] [-g [path]]

Databrowse: An Extensible Data Management Platform

optional arguments:

-h, --help	show this help message and exit
-s path, --setdataroot path	path to set new dataroot
-e, --openconfig	open cefdatabrowse config file
-g [path], --go [path]	open cefdatabrowse in a directory

On a fresh install a few configuration values will need to be changed. These include the dataroot directory and any third party software packages that need to be added to the Python path during operation. The CEFDatabrowse configuration file can be accessed using the command:

```
databrowse -e
```

This will open the configuration file in your default editor. From here the file will contain a default configuration like this:


```
[databrowse]
width = 800
height = 600
x = 0
y = 0
dataroot = C:/databrowse
limatix-qautils = C:/Users/nscheirer/Documents/dev/databrowse-utils/limatix-qautils
qautils = C:/Users/nscheirer/Documents/dev/databrowse-utils/QAutils
dataguzzlerlib = C:/Users/nscheirer/Documents/dev/dataguzzler-lib/python
```

From here you can set which directory you would like CEFDatabrowse to have access to. Keep in mind that CEFDatabrowse cannot access any files outside of this directory, therefore if you would like access to all files setting the dataroot to C:/ is acceptable. When a new dataroot is set it will be required that a directory called SOPs be created in the dataroot directory location. CEFDatabrowse will throw an error reminding you of this if you do not have a SOPs folder. The values following dataroot correspond to required third party packages that different plugins utilize. If you need that functionality replace these paths with the correct path to your installed copies. From this point CEFDatabrowse can be used via the defined usage above.

3 Getting Started

This section will provide some information on using the directory interface within Databrowse.

3.1 Features of Interest

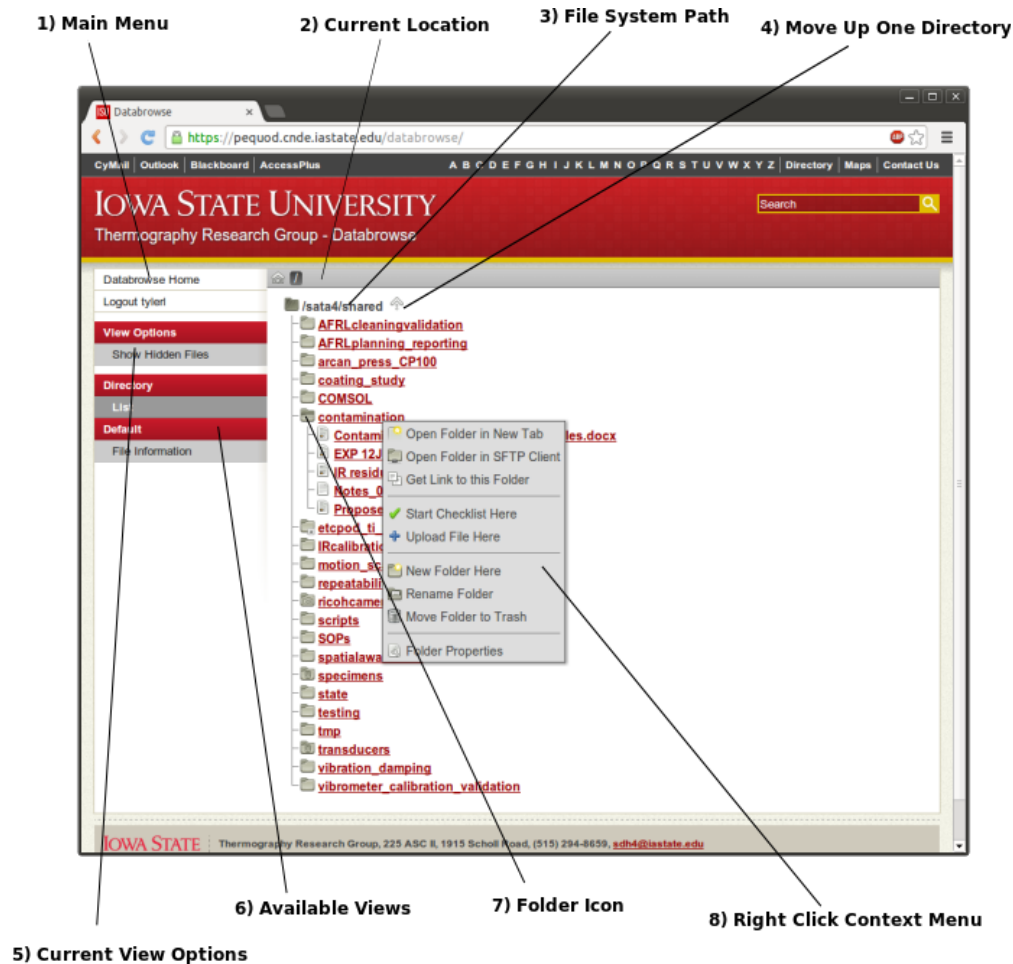

1. Main Menu

Contains links to jump to the data directory root, logout, and switch views

2. Current Location

Displays the location of the current view relative to the data directory root – click a folder name to jump to it

3. File System Path

Displays the absolute path of the current view in the file system

4. Move Up One Directory

Click this link to move up one directory until you reach the data directory root

5. Current View Options

Displays a list of options available to modify the current view

6. Available Views

Displays a list of various plugins that are capable of rendering the current file or folder and views available within those plugins

7. Folder Icon

Click this folder icon to expand the folder and dynamically load the contents of the sub folder

8. Right Click Context Menu

Right click the name of any folder or file in this view to display this context menu

3.2 Creating New Directories


The following steps can be used to create new directories from Databrowse:

1) Navigate to the Folder that will Contain the New Folder


Using the Databrowse interface, locate the folder that will contain the folder you wish to create.

2) Right Click Folder Name to Open Context Menu

Using the mouse, right click on the name of the folder you wish to contain your new folder. If the current view is the location you wish to create a new folder within, right click on the black text at the top. This is the shaded area shown on the figure to the right.

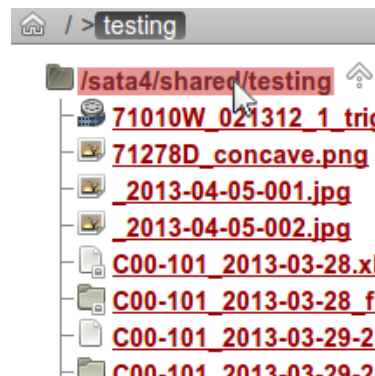


3) Select New Folder Here from the Context Menu

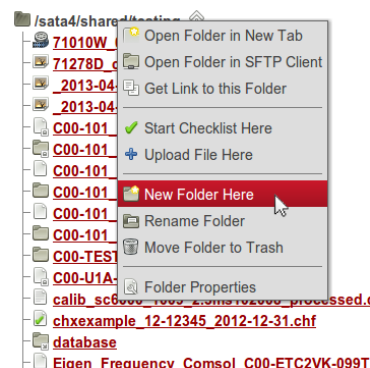

4) Type a Name and Click OK

Type a name for your new folder. Spaces are special characters that are permitted in folder names on your operating system may be used; however, their use is discouraged.

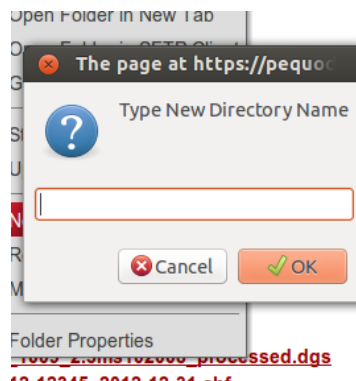

5) Finished!

You will be taken to the new folder automatically.

3.3 Uploading Files


The following steps can be used to upload files using Databrowse:

1) Navigate to the Folder in which you wish to Upload Files

Using the Databrowse interface, locate the folder that will contain the files you wish to upload.

2) Right Click Folder Name to Open Context Menu


Using the mouse, right click on the name of the folder in which you wish to upload files. If the current view is the location you wish to upload files to, right click on the black text at the top. This is the shaded area shown on the figure to the right.

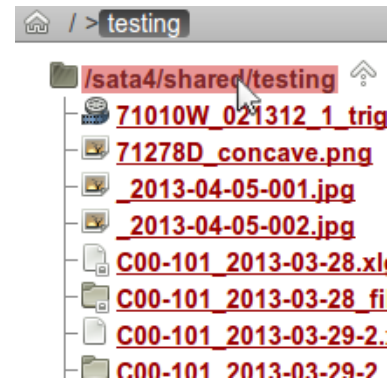

3) Select Upload Files Here from the Context Menu

4) Click Add Files to Open the File Selection Box

Inside the file selection box, you may hold Shift or Ctrl on the keyboard to select multiple files. Select Open once you have chosen the file(s) you wish to upload. Repeat for any additional files to be uploaded. The options you have selected will appear in the window.

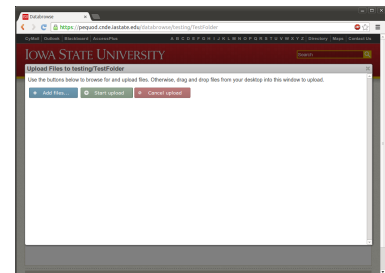


Alternatively, Drag and Drop Files into the Upload Window

Using Windows Explorer, Finder, or Nautilus, depending on your operating system, you may drag files into the open space on the upload window. The files will appear in the window.

5) Click Start Upload

You should see progress bars indicating the length of time remaining to upload your files. If uploading photos, you will see a preview thumbnail of the photo.

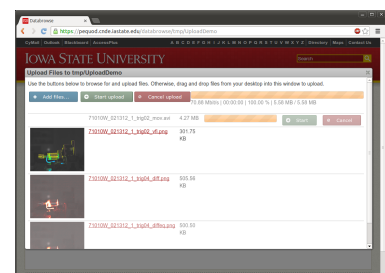

6) Click the X or Click Outside of the Upload Window to Close

The page will automatically refresh to show your new files. Do not close the window while uploading is still taking place or your upload will be interrupted.

4 Databrowse Configuration Options

This section provides information about the various configuration options that can be used to customize the low level functionality of Databrowse.

4.1 `databrowse_wsgi.conf`

The file `databrowse_wsgi.conf` (or a symbolic link to this file) must be contained inside the same folder as `databrowse.wsgi`. This file is a Python script that can be used to run code during the start of a web request to the server. It is ran within the context of the `web.support.py` file contained within `databrowse/support` in the source distribution. `web.support.py` is actually a Python class which is instantiated with the configuration from `databrowse_wsgi.conf` and the default values contained within the class at the start of a request. This code is ran directly during the call to `web.support.__init__`. Accordingly any Python code can be included within `databrowse_wsgi.conf` and it will be executed in this context.

The following options can be set in this file:

Site URL `self.siteurl` (*String*)

Default: `http://localhost/databrowse`

The full URL to Databrowse. This is used internally to build URLs. No trailing slash should be used if search engine optimized URLs are enabled. If you disable search engine optimized URLs, you should update this to list the full URL to `databrowse.wsgi`.

Resource URL `self.resurl` (*String*)

Default: `http://localhost/dbres`

The full URL to Databrowse static resource files. This is used internally to build URLs to static JavaScript, images, stylesheets, etc.

Logout URL `self.logouturl` (*String*)

Default: `http://localhost/logout`

The URL to use in order to trigger the web server's logout mechanism. This will need to be adjusted depending on your authentication method.

Icon Configuration File `self.icondbpath` (*String*)

Default: `os.path.join(os.path.dirname(databrowse.support.__file__), "iconmap.conf")`

The path to a file containing a ConfigParser-compatible configuration file associating file extensions with icons. See 4.2 for more information about the contents of this file.

Hidden File Configuration File `self.hiddenfiledbpath` (*String*)

Default: `os.path.join(os.path.dirname(databrowse.support.__file__), "hiddenfiles.conf")`

The path to a file containing a ConfigParser-compatible configuration file containing Glob syntax strings that can be used to hide certain files from view in the Databrowse directory views. See 4.3 for more information about the contents of this file.

Search Engine Optimization `self.seo_urls` (*Boolean*)

Default: `True`

When set to true, Databrowse will write out URLs that have been search engine optimized. Instead of the URL `http://localhost/databrowse.wsgi?path=/SOPs`, the URL would be written as `http://localhost/databrowse/SOPs`. Usage of this feature requires `mod_rewrite` to be enabled with the options as suggested in the Installation instructions. It is strongly encouraged to keep this feature turned on.

Other options can also be set in this file using the same format, if desired, for use in Databrowse plugins.

4.2 iconmap.conf

This file generally does not need to be modified, unless you wish to build or modify plugins. This file tracks the mapping between file types and icons. These icons are used to provide a graphical representation of a file type within Databrowse. The file uses a Python ConfigParser module compatible format. There are two sections that can be used inside this file: `[Content-Type]` and `[Extension]`. Content types are defined as specified in RFC 2045 and are determined utilizing libmagic. File extensions are matched against the portion of the filename after the final period. The content-type or extension becomes the parameter name in `iconmap.conf` and the parameter value is the name of the file containing the icon served at the URL `icons/` relative to the Databrowse resources URL (e.g. `http://localhost/dbres/icons/folder.png`).

The following is an example of how this file is laid out:

```
[Content-Type]
inode/directory=folder.png
application/x-directory=folder.png

[Extension]
txt=text-x-generic.png
html=text-html.png
```

All files that are unable to be matched with a configuration option are given the icon `unknown.png`.

4.3 hiddenfiles.conf

This file can be used to hide files based on their file name. Glob expressions are permitted. This file uses a Python ConfigParser module compatible format. There are two sections that can be used inside this file: `[Hidden]` and `[Shown]`. The hidden section is used to define file names or glob expressions for filenames that should be hidden from the Databrowse directory view. The shown section can be used to override glob expressions contained in the hidden section for explicit file names. Glob expressions are also permitted here. The parameter name is used only for comment and will be ignored by Databrowse. The parameter value will be used by Databrowse to filter file lists in the directory view.

The following is an example of how this file is laid out:

```
[Hidden]
backup_files=*~
hidden_files=.*
generic_backup_files=*.bak

[Shown]
my_special_backup_file=MyBackup*.bak
```

5 Databrowse Library

The core components of Databrowse that are most directly responsible for building XML representations of files and directories are packaged into Python modules. An interface has been provided that enables Python scripts to leverage the power of Databrowse to quickly obtain XML representations for any file. This is particularly useful in the context of processing scripts written in Python. This is also especially useful in contexts in which the Databrowse plugin responsible for handling a particular type of file provides additional information or meta data that wouldn't otherwise be easily accessible working with the file directly.

The Databrowse library can be imported in Python with the following code:

```
from databrowse.lib import db_lib as dbl
```

The library contains one function named `GetXML`. This function will return the XML representation of a file, as produced by Databrowse.

```
function dbl.GetXML(filename, output=dbl.OUTPUT_ELEMENT, **params)
```

Aguments

filename String containing a relative or absolute path to the file of interest

output Determines the type of output to be returned from the function

dbl.OUTPUT_ELEMENT returns an LXML etree.Element

dbl.OUTPUT_ETREE returns an LXML etree.ElementTree

dbl.OUTPUT_STRING returns a string containing the XML

dbl.OUTPUT_STDOUT prints the XML string to stdout and returns nothing

****params** A variable number of optional parameters that are treated the same way as query string values that would be POST or GET to the web server when Databrowse is being used from the web. Used to pass in various options into plugins.

Usage

```
>>> from databrowse.lib import db_lib as dbl
>>> dbl.GetXML('/tmp/emptyfile', output=dbl.OUTPUT_STDOUT)
<default:default>
  <filename>emptyfile</filename>
  <path>/tmp</path>
  <size>0.0 byte</size>
  <mtime>Tue Sep  3 10:12:40 2013</mtime>
  <ctime>Tue Sep  3 10:12:40 2013</ctime>
  <atime>Tue Sep  3 10:12:42 2013</atime>
  <contenttype>text/plain</contenttype>
  <permissions>-rw-rw-r--</permissions>
  <owner>user:user</owner>
</default:default>
```

It is also worth mentioning the presence of the function `DebugGetXML` which operates identically to `GetXML`, however, it launches a Python debugger session, enabling the user to step through the code line by line.

6 Databrowse Plugin Format

This section provides documentation on nature of Databrowse plugins. It is intended to provide enough information that skilled users may be able to build or customize their own Databrowse plugins, but it is not intended to be a tutorial. Such material will be produced for Databrowse v1.0.

6.1 Plugin File Structure

Individual plugins are Python packages contained within the `databrowse.plugins` namespace. All plugin names should start with the prefix `db_` (this is a legacy requirement and will be removed in Databrowse v1.0). The file structure of the package should be as follows:

```
db_plugin_name\  
    __init.py__  
    db_plugin_name.py  
    dbs_stylesheet_one.xml  
    dbs_stylesheet_two.xml  
    handlers.py
```

6.2 File Contents

This section will detail the required contents for each file in the plugin. For the purposes of providing an example, let's say we wish to construct a plugin that simply displays the name of a file to the web page. This section will display the necessary file contents needed to produce such a plugin.

6.2.1 `__init.py__`

This file does not require any contents. Its presence is a trigger to the Python interpreter to search for modules. However, you may find it useful to place a docstring inside this file for usage from the Python console and for documentation purposes. The presence of a copyright statement is strongly encouraged. This file can also be used to initialize variables, though, this usage is strongly discouraged in the context of the WSGI server.

6.2.2 `db_plugin_name.py`

This file must be named with the same name as the folder name. The top of the file should contain a copyright statement. The usage of a docstring at the top of the file is also recommended for identification purposes.

This file must contain a class derived from `databrowse.support.renderer_support.renderer_class` and named identically to the file name. It must have several class variables defined:

`_namespace_uri` A string containing the fully qualified URL referring to the XML namespace for the plugin. Existing Databrowse plugins follow the format `http://thermal.cnde.iastate.edu/databrowse/plugin_name`.

`_namespace_local` A string containing the local abbreviated form of the namespace. You should take care to avoid conflicts with other namespace prefixes.

Several other class variables are automatically initialized to defaults, but may be overridden by declaring them as class variables to your inherited class:

`_default_content_mode` A string identifying the default content mode to be loaded when no content mode is specified. The content mode `full` is frequently used. This serves more as an internal identifier for use within the plugin; however, the content mode `raw` has special meaning. This content mode will prevent the main WSGI script from outputting any content to the user. As a result, you can use the `raw` content mode to serve binary content to the web browser. You must manually set the output headers and return content as appropriate from your plugin per PEP 3333.

_default_style_mode A string identifying the name of the stylesheet that should be applied to the content output by default if one is not specified by the end user. This string should not include the **db_** prefix. See Section 6.2.3 for more information.

_default_recursion_depth A number indicating how deep a plugin operating on a directory should recurse down a directory tree by default. This value is not used by plugins that do not operate on directories but should still be set.

An additional set of class variables will be initialized to contain references to objects used by other portions of Databrowse that may be of convenience or use here:

_relpath A string initialized to the file path of the currently requested file relative to the data root with a leading forward slash. These paths are used by Databrowse internally to represent full file paths inside URLs.

_fullpath A string initialized to the absolute file path of the currently request file. This path should be used to access the file.

_web_support A reference to the instantiated **web_support** class, which contains information about the web server request and the current configuration of Databrowse.

_handler_support A reference to an instantiated **handler_support** class, which provides functionality for determining which types of plugins are capable of operating on a particular file. It also provides support functionality to determine the icon used to represent a file.

_caller A string that identifies the plugin responsible for the call to this plugin. When set to **databrowse**, the call to the plugin is the result of a direct request from the Databrowse WSGI application. Otherwise, this string will contain the name of the plugin making the recursive call to get content from the plugin. This is most frequently the Directory plugin, which can query a plugin for information about a file. This information can then be displayed by the Directory plugin as a preview of a file's contents.

_handlers A list of other plugins that are capable of working with this particular type of file, in order of precedence.

The class must contain at least one function – **getContent**. With the exception of the scenario in which the content mode is set to **raw**, this function must return via with an LXML etree.Element, **None**, or raise an exception. The keyword **None** should only be returned when the plugin is not called by Databrowse, but rather is being called by another plugin. This is normally the desired behavior, unless your plugin needs to return content to a directory plugin. The example contents of **db_plugin_name.py** displayed on the following page will produce an XML document containing the filename of the requested file. The format of the generated XML document is

```
<pn:pn xmlns:pn="http://thermal.cnde.iastate.edu/databrowse/plugin_name">
    <filename>some_file_name.txt</filename>
</pn:pn>
```

if the plugin were called with the filename **some_file_name.txt**. This example plugin is trivial, but provides the framework needed to produce XML representations of any type of file.

```

#!/usr/bin/env python
#####
## Databrowse: An Extensible Data Management Platform ##
## Copyright (C) 2012-2015 Iowa State University ##
## ##
## This program is free software: you can redistribute it and/or modify ##
## it under the terms of the GNU General Public License as published by ##
## the Free Software Foundation, either version 3 of the License, or ##
## (at your option) any later version. ##
## ##
## This program is distributed in the hope that it will be useful, ##
## but WITHOUT ANY WARRANTY; without even the implied warranty of ##
## MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the ##
## GNU General Public License for more details. ##
## ##
## You should have received a copy of the GNU General Public License ##
## along with this program. If not, see <http://www.gnu.org/licenses/>. ##
#####
""" plugins/db_plugin_name/db_plugin_name.py - The Plugin Name Plugin """

from lxml import etree
from databrowse.support.renderer_support import renderer_class

class db_plugin_name(renderer_class):
    """ The Plugin Name Plugin - This plugin is an example """

    _namespace_uri = "http://thermal.cnde.iastate.edu/databrowse/plugin_name"
    _namespace_local = "pn"
    _default_content_mode = "full"
    _default_style_mode = "default_view"
    _default_recursion_depth = 2

    def getContent(self):
        if self._caller != "databrowse":
            return None
        else:
            if self._content_mode == "full":
                xmlroot = etree.Element('{%s}%s' % (self._namespace_uri, self._namespace_local),
                                         nsmap=self.nsmap)
                xmlchild = etree.SubElement(xmlroot, "filename", nsmap=self.nsmap)
                xmlchild.text = os.path.basename(self._fullpath)
                return xmlroot
            else:
                raise self.RendererException("Invalid_Content_Mode")
        pass
    pass

```

6.2.3 dbs_stylesheet_name.xml

Plugins can contain any number of XSLT stylesheets that are used to transform the XML content produced by the plugin for a particular file into HTML for display on the web browser. These stylesheet files must be named with the prefix `dbs_` and end with the extension `.xml`. Databrowse will automatically locate all stylesheets located within the plugin named in this format and display them as options to the user in the Databrowse menu. Databrowse will also search for stylesheets inside the data folders placed inside of the `.databrowse/stylesheets/db_plugin_name` folder, enabling the end user to write stylesheets for custom views relevant to a particular set of data.

This file does not contain a complete XSLT stylesheet. Rather, it contains XSLT template snippets which will be combined with other snippets from other plugins that may be producing content displayed on the page at any given time. This is particularly important in the context of a directory plugin that is displaying representations of many different files contained within that directory.

This file also does not produce a complete web page. Rather, it should be written such that it writes out HTML snippets, which will be placed in the appropriate location on the web page.

The construction of an XSLT template and the usage of HTML is outside of the scope of this document. Please refer to one of many tutorials on the topics available on line.

The following is an example of an XSLT stylesheet that could be used to transform the content produced by the discussed plugin into HTML.

```
<xsl:template xmlns="http://www.w3.org/1999/xhtml"
  xmlns:pn="http://thermal.cnde.iastate.edu/databrowse/plugin_name"
  xmlns:xsl="http://www.w3.org/1999/XSL/Transform" match="pn:pn" mode="full">
  <h1><xsl:value-of select="filename" /></h1>
</xsl:template>
```

6.2.4 handlers.py

Plugins must contain a `handlers.py` file. This file should contain one function with a name starting with the prefix `dbh_`. The name of this function will be used to determine the order of precedence in ascending order in the case in which multiple plugins can operate on a file. This function receives three parameters:

path A string containing the full path to the file that has been requested.

contenttype A string containing the RFC 2045 compliant MIME type.

extension A string the portion of the filename after the final period or empty if there is no period in the filename.

The return from this function should either be the name of a plugin capable of responding to a request for the given file or `False` otherwise.

Returning to the previous example, a `handlers.py` file for our example plugin might look like the file displayed below. The displayed code will result in this plugin being able to respond to any file with the `.txt` extension.

```
#!/usr/bin/env python
#####
## Databrowse:  An Extensible Data Management Platform      ##
## Copyright (C) 2012-2015 Iowa State University            ##
##                                                         ##
## This program is free software: you can redistribute it and/or modify ##
## it under the terms of the GNU General Public License as published by ##
## the Free Software Foundation, either version 3 of the License, or    ##
## (at your option) any later version.                            ##
##                                                         ##
## This program is distributed in the hope that it will be useful,      ##
## but WITHOUT ANY WARRANTY; without even the implied warranty of      ##
## MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE.  See the      ##
## GNU General Public License for more details.                  ##
##                                                         ##
## You should have received a copy of the GNU General Public License    ##
## along with this program.  If not, see <http://www.gnu.org/licenses/>. ##
#####
""" plugins/handlers/dbh_plugin_name.py - Handler for the Plugin Name plugin """

def dbh_plugin_name(path, contenttype, extension):
    """ Plugin Name Handler - Responds to *.txt files"""
    if extension == "txt":
        return "db_plugin_name"
    else:
        return False
```


7 Included Databrowse Plugins

This section details the various plugins that are included with Databrowse for working with a variety of different file types. Documentation in this section is a work in progress. For the time being, simple descriptions of each plugin have been provided.

7.1 Checklist Editor

The checklist editor plugin is one of several Databrowse plugins that were constructed to create and modify data. Checklists (*.chx for checklist templates and *.chf for filled checklists) are XML files that contain the necessary information to document laboratory procedures. The complete details of the file format are outside of the scope of this document. However, this plugin enables the user to open a checklist template, fill out a checklist, and save the filled checklist to the file system.

An experimental tool for modifying a checklist is being developed as well. This tool utilizes the Axel XML JavaScript library (<http://ssire.github.io/axel/>) to provide an interface for editing the file in the web browser.

This plugin is not intended to be used from within the Databrowse library interface.

The screenshot shows the Databrowse web interface with the 'Checklist Editor' plugin active. The left sidebar contains navigation links: 'Checklist Editor', 'File Checklist', 'Generate XML File', 'XML Preview', 'Plain Text File', 'Edit Text File', 'Plain Text Preview', and 'Default'. The main content area displays a checklist titled 'Repeatability Startup'. The checklist has columns for 'Checklist', 'Specimen', 'Performed by', 'Date', and 'Destination'. The 'Checklist' column contains seven items, each with a description and a 'Last Updated' timestamp. The 'Specimen' column contains 'C801-TEST'. The 'Performed by' column contains 'Tyler'. The 'Date' column contains '2015-08-28'. The 'Destination' column contains 'trng'. The checklist items are: 1. Open the Run Sequence Document or Obtain Printout, 2. Determine the Run Number to be Completed, 3. Ensure the Performed By Field is Filled Above, 4. Prepare the Thermal Camera and Vibrometer, 5. Ensure Dataguzzler is not already Running, 6. Turn on the Waveform Generator, Laser Vibrometer, and Thermal Camera, and 7. Start Dataguzzler.

7.2 Checklist Viewer

The checklist viewer plugin displays filled checklist files (*.chf). Checklists (*.chx for checklist templates and *.chf for filled checklists) are XML files that contain the necessary information to document laboratory procedures. The complete details of the file format are outside of the scope of this document. However, this plugin will display the status of all checklist items, along with displaying time stamp information and notes.

The screenshot shows the Databrowse web interface with the 'Checklist Viewer' plugin active. The left sidebar contains navigation links: 'Checklist Viewer', 'View Filled Checklist', 'Checklist Editor', 'Edit Checklist', 'Fill Out Checklist', 'Generate XML File', 'XML Preview', 'Plain Text File', 'Edit Text File', 'Plain Text Preview', and 'Default'. The main content area displays a checklist titled 'Repeatability Shutdown'. The checklist has columns for 'Checklist', 'Specimen', 'Performed by', 'Date', and 'Destination'. The 'Checklist' column contains nine items, each with a description and a 'Last Updated' timestamp. The 'Specimen' column contains 'C801-TEST'. The 'Performed by' column contains 'Tyler'. The 'Date' column contains '2015-08-28'. The 'Destination' column contains 'trng'. The checklist items are: 1. Move Z Stage (if necessary), 2. Turn Pressure Off, 3. Turn off the IR Camera, 4. Turn off the Laser Vibrometer, 5. Replace Lens Caps on IR Camera and Laser Vibrometer, 6. Turn off the Waveform Generator, 7. Turn off the Power Amplifier, 8. Quit Dataguzzler, and 9. Press One Red Emergency Shutdown Button.

7.3 Datacollect v1 Viewer

Datacollect is a tool that enables the automation of data collection processes. Version 1 of Datacollect provides an interface for the entry of experimental parameters and an interface to aid in automating use of data acquisition tools. This plugin is responsible for displaying the data saved by Datacollect in a meaningful fashion. It can display data in both a log-style format and in a tabular form, capable of being sorted, filtered, and searched.

Database Home
Logout Titled

Datacollect V1 Viewer
Tabular Summary
Log View
Generic XML File
XML Preview
Plain Text File
Grid Test File
Plain Text Preview
Default
File Information

Run Summary

Run Index: 1

Date: 10/25/07

Baseline:

Description:

Variable: amplitude

Label: etcpod

Trigger Data

Show 10 entries

Air Cylinder	Amplitude	Clamp Offset	Clamp Torque	Coating	Couplant	Crat
McMaster 1.5in	4.0000000	10.5	28	icepaint	ETC standard paper	0.02
McMaster 1.5in	3.0000000	10.5	28	icepaint	ETC standard paper	0.03
McMaster 1.5in	3.0000000	10.5	28	icepaint	ETC standard paper	0.09
McMaster 1.5in	1.5000000	10.5	28	icepaint	ETC standard paper	0.02
McMaster 1.5in	2.2000000	10.5	28	icepaint	ETC standard paper	0.01
McMaster 1.5in	1.5000000	10.5	28	icepaint	ETC standard paper	0.00
McMaster 1.5in	2.2000000	10.5	28	icepaint	ETC standard paper	0.02

1 (1 entries shown)

Showing 1 to 7 of 7 entries

ShowHide Fields

Air Cylinder
Coating
Droaine Stress
Hil
Resonant Frequency
TimeJamp
Transducer Serial
Amplitude
Couplant
Electrical Energy
Measurement Type
Run Basename
Transducer Force
Transducer Temperature
Clamp Offset
Crat Heating
Excitation
Notes
Set File
Transducer Offset
Trigger Count
Clamp Torque
DGS File
HSE
FilesInCh
Taken By
Transducer Pressure

7.4 Datacollect v2 Viewer

Datacollect is a tool that enables the automation of data collection processes. Version 2 of Datacollect provides an interface for the entry of experimental parameters, an interface to aid in automating use of data acquisition tools, and an interface for managing experimental processes and procedures. This plugin is responsible for displaying the data saved by Datacollect in a meaningful fashion. It can display data in both a log-style format and in a tabular form, capable of being sorted, filtered, and searched.

- Database Home
- Logout
- Log View
- Raw Data
- Generic XML File
- XML Preview
- Plain Text File
- Edit Text File
- Plain Text Preview
- Details
- File Information

[Go](#)
[>](#)
[UT0vtrb035](#)
[>](#)
[empirical_modeling](#)
[>](#)
[UGdata](#)
[>](#)
[C14-UTCA-002X_turboisy_2014-09-18.xml](#)

Datacollect Experiment Log

Experiment Summary

Item	Details
Specimen:	C14-UTCA-002X
Performed By:	Leshaighe
Date:	2014-09-18:2014-09-19:2014-10-11
Data Folder:	C14-UTCA-002X_turboisy_2014-09-18_files
Goal:	
Experiment Notes:	

Experiment Photographs

Configuration

[/database/SOPs/OpticalScanning/support/crackopening.docx](#)

```

# Note: datacollect config files are passed through an m4 preprocessor
# so you can use e.g. "m4_include([[includefile.doc]])" to include
# another file. As an alternative you can manually load multiple config
# files before opening your experiment log, either from the
# datacollect command line or from the configuration gui

# Call dglo.starttime() ONLY if this config needs to talk to datagazer!
dglo.starttime() # multiple calls to dglo.starttime() are OK

# call createparameter() if #bus (i.e. python/matlab/etc.) access to the parameter
# database is likely to be needed.
createparameter() # Multiple calls to createparameter() are OK
  
```

Measurement Summary

Item	Value
Measurement #	0
Specimen	C14-UTCA-002X
Date	2014-09-18
Performed By	Vadsi
Checklist Name	vibroneorientation_press_V1_2_2014-05-16
Checklist Title	Vibrotomography New Specimen or Orientation testing in the press
Checklist Link	C14-UTCA-002X_turboisy_2014-09-18_C14-UTCA-002X_neworientation_press-0009.cdf
DGS File	C14-UTCA-002X_turboisy_2014-09-18_C14-UTCA-002X_neworientation_press-

7.5 Dataguzzler Data File

Dataguzzler is an open source extensible data acquisition platform – providing a high speed mechanism for data capture, storage, and visualization. Dataguzzler also provides bindings enabling simple integration with data collection automation scripts and other related programs. The Dataguzzler data file plugin for Databrowse enables the user to view the data contained within Dataguzzler binary data files (*.dgs, *.dgd, *.dga, *.dgz). The plugin utilizes Matplotlib for Python to produce realtime visualizations of data, in addition to displaying all of the associated meta data. The plugin can also export data from Dataguzzler data files into several other common file formats, including CSV and MAT. The plugin is also capable of exporting videos from the appropriate types of data.

Database Home
Logout taylor

C14-040-004g_FlawB_tortuosity_2015-02-17_files > C14-040-004g_FlawB_tortuosity_2015-02-17_C14-040-004g_testfreq-0018.dgs

Dataquizzer Data File
Examine File Contents
View Raw Structure
Generic Binary File
Binary Preview
Default
File Information

C14-040-004g_FlawB_tortuosity_2015-02-17_C14-040-004g_testfreq-0018.dgs

Dataquizzer Snapshot

Waveform Count: 31
Time Stamp: 2015-02-20T04:01:32-0500
Scalar: Energy HAE
Waveforms: 434.87 Joules 0.106095 meters²/second


Waveforms of Interest

View Waveform:

Azb	ArtGen	DirStack	DynamicStrain
Energy	Force	HAE	Rmg
IStack	Itt	Pressure	Vib
VibAct	VibQm	VibCompensate	VibAct
VibSensep	VibSensepFFT	Vibm	VibmFds
VibmFitting	Vibrometer	Vibrometer2	VibmCompensate
VibmACVoltage	VibmCurrent	VibmDebias	VibmVoltage
IR	hertzian	cath	

7.6 Dataguzzler Settings File

Dataguzzler is an open source extensible data acquisition platform – providing a high speed mechanism for data capture, storage, and visualization. Dataguzzler also provides bindings enabling simple integration with data collection automation scripts and other related programs. The Dataguzzler set file plugin for Databrowse enables the user to examine the contents of Dataguzzler settings files (*.set).

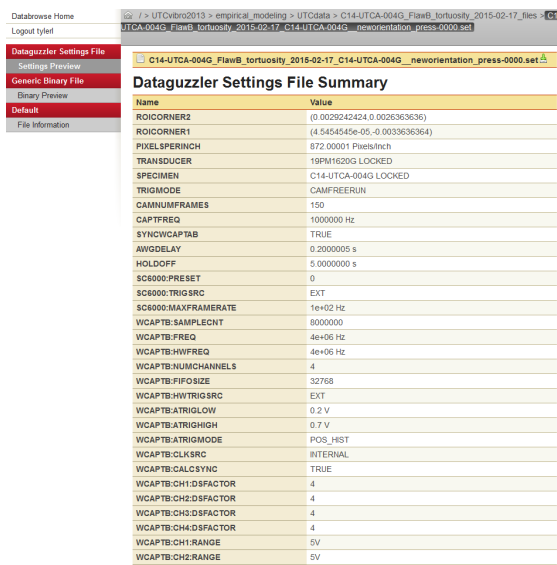

The screenshot shows the Databrowse application window. The 'Database Home' is set to 'C:\> \UTCORR013 > empirical_modeling > UTCData > C14-UTCA-004G_FlawB_tortuosity_2015-02-17_files'. The 'Logon User' is 'UTC004G_FlawB_tortuosity_2015-02-17\C14-UTCA-004G_neworientation_press-0000.sci'. The 'Dataguzzler Settings File' is selected in the left sidebar. The main pane displays the 'Dataguzzler Settings File Summary' table.

Name	Value
ROICORNER2	(0.0029242424,0.0026363636)
ROICORNER1	(4.5454545e-05,-0.0036363634)
PIXELSPERINCH	872.00001 Pixels/inch
TRANSDUCER	19PM1620G LOCKED
SPECIMEN	C14-UTCA-004G LOCKED
TRIGMODE	CAMFREERUN
CAMNUMFRAMES	150
CAPTFRQ	1000000 Hz
SYNCGCAPTAB	TRUE
AWGDELAY	0.2000005 s
HOLDOFF	5.0000000 s
SC6000:PRESET	0
SC6000:TRIGSRC	EXT
SC6000:MAXFRAMERATE	1e+02 Hz
VCAPTIB:SAMPLECHIT	8000000
VCAPTIB:FREQ	4e+06 Hz
VCAPTIB:NUMCHANNELS	4
VCAPTIB:FIPOSIZE	32768
VCAPTIB:HWTRIGSRC	EXT
VCAPTIB:ATRIGLOW	0.2 V
VCAPTIB:ATRIGHIGH	0.7 V
VCAPTIB:ATRIGMODE	POS_HST
VCAPTIB:CLKSRC	INTERNAL
VCAPTIB:CALCSYNC	TRUE
VCAPTIB:CH1:DSFACTOR	4
VCAPTIB:CH2:DSFACTOR	4
VCAPTIB:CH3:DSFACTOR	4
VCAPTIB:CH4:DSFACTOR	4
VCAPTIB:CH1:RANGE	5V
VCAPTIB:CH2:RANGE	5V

7.7 Data Table

The data table plugin is a tool designed to aid in the collection and compilation of large data sets from multiple sources. This plugin operates on *.tbl files. The exact schema of the *.tbl file is outside of the scope of this document; however, the tbl file is an XML file which instructs the Data Table plugin where to find data files and how to select data from them. The end result is a new XML document that combines all of this data together. This XML document can then be rendered by Databrowse in the web browser in the form of a searchable, sortable, and filterable table. This table can be exported to CSV from the web browser as well. Additionally, this plugin can prove to be very powerful when utilized inside of a processing script using the Databrowse library, enabling rapid development of data processing scripts that can operate on large sets of data in real time.

Database Home	C:\> \UTCORR013 > empirical_modeling > UTCData > C14-UTCA-004G_FlawB_tortuosity_2015-02-17_files					UTCORR013 > empirical_modeling > UTCData > C14-UTCA-004G_FlawB_tortuosity_2015-02-17_files	
Logon User	UTC004G_FlawB_tortuosity_2015-02-17\C14-UTCA-004G_neworientation_press-0000.sci					UTC004G_FlawB_tortuosity_2015-02-17\C14-UTCA-004G_neworientation_press-0000.sci	
Processed Data							
Database	File Name	Frequency	Default	Actual	Units	Source	Units
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-01.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-02.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-03.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-04.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-05.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-06.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-07.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-08.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-09.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-10.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-11.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-12.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-13.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-14.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-15.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-16.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-17.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-18.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-19.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-20.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-21.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-22.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-23.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-24.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-25.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-26.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-27.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-28.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-29.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-30.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-31.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-32.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-33.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-34.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-35.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-36.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-37.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-38.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-39.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-40.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-41.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-42.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-43.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-44.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-45.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-46.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-47.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-48.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-49.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-50.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-51.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-52.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-53.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-54.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-55.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-56.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-57.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-58.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-59.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-60.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-61.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-62.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-63.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-64.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-65.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-66.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-67.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-68.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-69.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-70.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-71.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-72.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-73.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-74.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-75.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-76.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-77.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-78.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-79.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-80.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-81.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-82.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-83.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-84.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-85.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-86.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-87.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-88.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-89.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-90.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-91.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-92.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-93.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-94.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-95.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-96.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-97.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-98.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-99.tbl	1000000	1000000	1000000	Hz	1000000	Hz
C14-UTCA-004G	UTC004G_FlawB_tortuosity_2015-02-17-100.tbl	1000000	1000000	1000000	Hz	1000000	Hz

7.9 Directory



The generic directory plugin is capable of scanning a folder and building an XML representation of the contents of that folder. It is also capable of recursively building a representation of entire directory trees. In the web interface, this is displayed as the default file browser view. The web browser interface has a number of useful features that enable the user to perform operations on files. See Section 3 for additional information about this interface.

7.10 File Operations

The file operations plugin is a special plugin that is normally not visible from the web or library interface. This plugin enables other plugins to perform special operations on files, particularly those that do not yet exist.

7.11 Generic Binary File

The generic binary file plugin is capable of reading out data in a hex and ASCII format from a binary data file. In the web interface, generic information is displayed about the file, along with an interactive hex viewer tool. The viewer will only pull in a small chunk of the file at a time, since binary files can be rather large. AJAX requests are made by the web browser enabling the user to browse through the file interactively in real time.

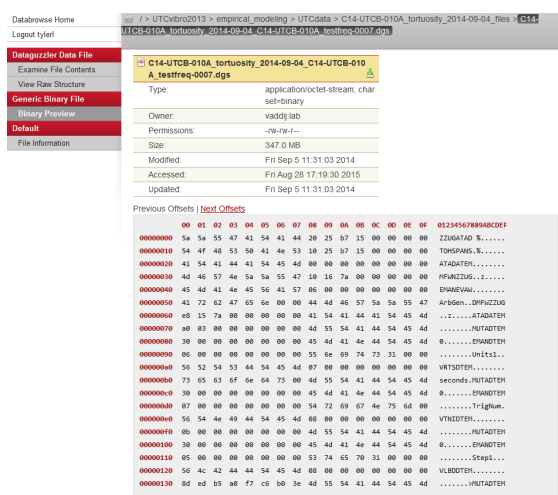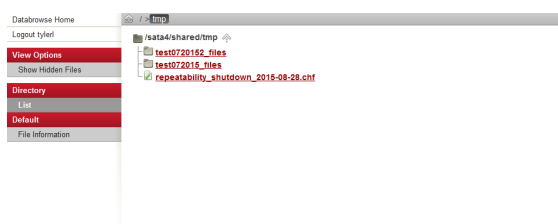

7.12 Generic HDF5 File

The HDF5 viewer plugin is a work-in-progress plugin that will simply display the internal structure of an HDF5 file at this time.

The screenshot shows the Databrowse application with the 'Generic HDF5 File' plugin selected. The file 'cubicSingleQuaxed.dream3d' is displayed with its metadata. The main content area shows the internal structure of the HDF5 file as XML, including details about the schema, dimensions, and data types.

7.13 Generic WSGI Application

The generic WSGI application plugin will enable a WSGI compliant Python script to be ran inside the context of Databrowse. This is very useful for scenarios in which an end user wishes to quickly create an interactive web application within the context of a particular set of data or for a task-oriented purpose without the need to create an infrastructure in which that application will run.


7.14 Generic XML File

The generic XML file viewer will display a textual view of an XML file, along with basic file information.

The screenshot shows the Databrowse application with the 'Generic XML File' plugin selected. The file 'etcpod_05-400_102307_1.xml' is displayed with its metadata. The main content area shows the XML content of the file, which appears to be a configuration or data file for a system.

7.15 Image Viewer


The image viewer plugin will display an image, in addition to all internal metadata and EXIF tags associated with an image. The plugin also supports the ability to resize images in real time for use by other plugins as well. It will operate on any file type supported by the Python Imaging Library, including *.png, *.jpg, *.bmp, and many more.

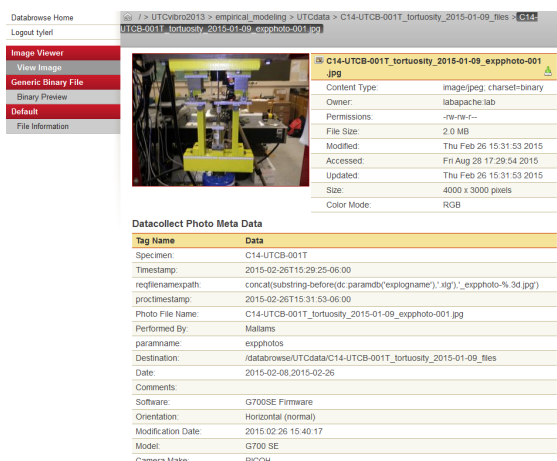

Tag Name	Data
Specimen	C14-UTCB-001T
Timestamp	2015-02-26T15:29:25-06:00
refilenamespath	concat(substring-before(dc:parandb('explogname'),'.xlg'),'_expphoto-%.3d.jpg')
proctimestamp	2015-02-26T15:31:53-06:00
Photo File Name	C14-UTCB-001T_tortuosity_2015-01-09_expphoto-001.jpg
Performed By	Mallams
paramname	expphotos
Destination	/databrowse/UTdata/C14-UTCB-001T_tortuosity_2015-01-09_files
Date	2015-02-08,2015-02-26
Comments	
Software	G700SE Firmware
Orientation	Horizontal (normal)
Modification Date	2015-02-26 15:40:17
Model	G700 SE
Camera Make	RICOH

7.16 Mercurial Repository


Mercurial is a version tracking tool. The Mercurial Repository Databrowse plugin is a work-in-progress that will display warning information to users about uncommitted changes and other potential concerns associated with a Mercurial repository. This plugin overrides the default Directory plugin, though, all other functionality available from the directory plugin is still available in this context.

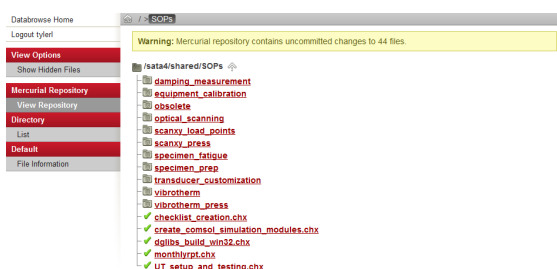

- data4/shared/SOPs
 - damping_measurement
 - equipment_calibration
 - obsolete
 - optical_scanning
 - scanxy_load_points
 - scanxy_press
 - specimen_fatigue
 - specimen_press
 - transducer_customization
 - vibrotherm
 - vibrotherm_press
 - ✓ checklist_creation.chx
 - ✓ create_console_simulation_modules.chx
 - ✓ dglibs_build_win32.chx
 - ✓ monthlrpt.chx
 - ✓ UT_setup_and_testing.chx

7.17 Movie Viewer


The movie viewer plugin is an experimental plugin that will attempt to stream video content to the users web browser. This requires the usage of the appropriate plugins and codecs on the end user computer. This plugin is also capable of extracting preview frames from videos. It will operate on any video format supported by FFMPEG.

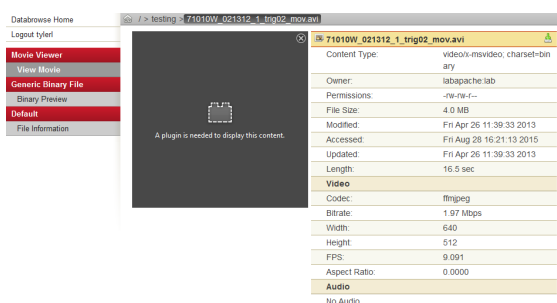

Tag Name	Data
Content Type	video/mpeg, charset=bin
Owner	labapache:lab
Permissions	-rwxr-xr-x
File Size	4.0 MB
Modified	Fri Apr 26 11:39:33 2013
Accessed	Fri Aug 26 16:21:13 2015
Updated	Fri Apr 26 11:39:33 2013
Length	16.5 sec
Video	
Codec	ffmpeg
Bitrate	1.97 Mbps
Width	640
Height	512
FPS	9.091
Aspect Ratio	0.0000
Audio	No Audio


7.18 Multimedia Directory

The multimedia directory plugin will display thumbnail previews of image, video, and other multimedia content. It will override the default directory plugin when the majority of the content inside of a folder can be displayed with a thumbnail preview. Clicking on a thumbnail in the multimedia directory plugin will display a larger version of the image without leaving the page.

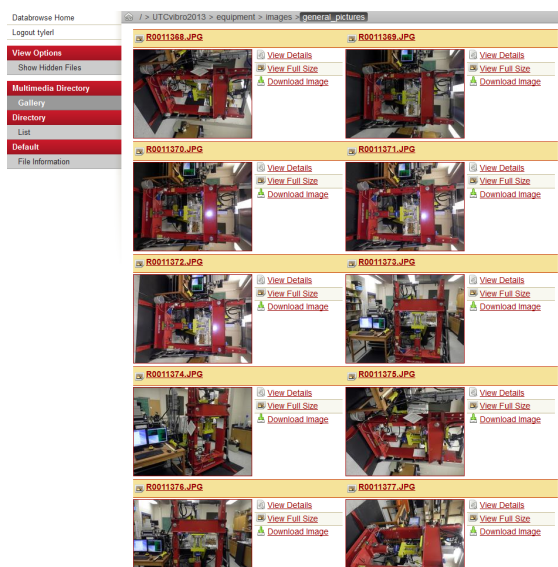


7.19 Office Viewer

The office viewer plugin provides an in-browser PDF preview of word processor documents, spreadsheets, presentations, and other document files. This plugin operates on any file format supported by LibreOffice, including *.doc, *.ppt, *.xls, *.odt, *.odp, *.ods, and many others.

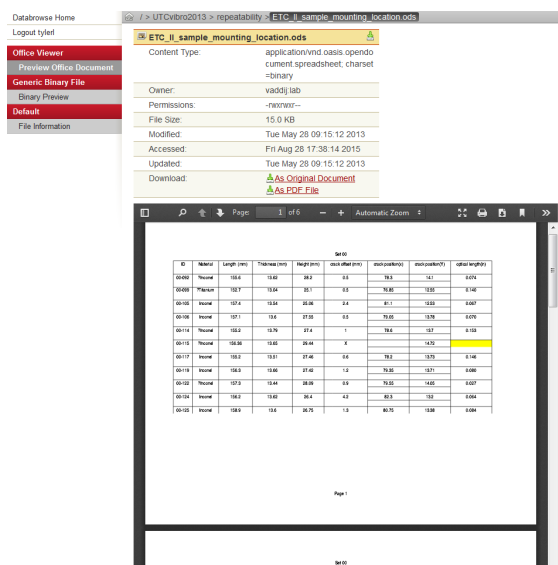


7.20 PDF Viewer

The PDF viewer plugin will display a PDF file for the user inside the web browser.

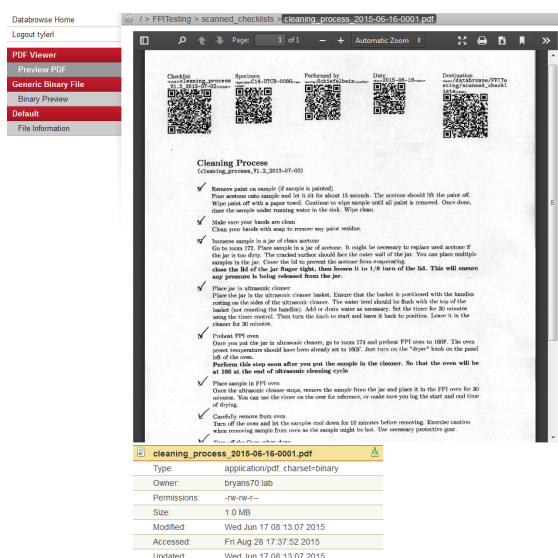


7.21 Plain Text File

The plain text file plugin displays the contents of any plain text file inside the web browser. It also provides a convenient interface for editing plain text files right in the web browser. Backups will automatically be saved by this plugin during a save operation. It will operate on all files with the content-type of plain/text.

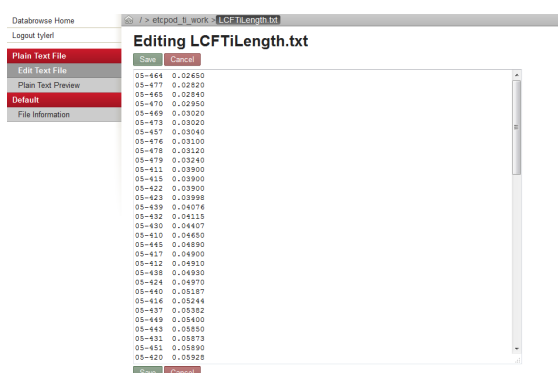

7.22 SolidWorks Viewer

The SolidWorks viewer plugin provides a quick low resolution image preview of a SolidWorks file in the web browser. This plugin is also capable of providing a thumbnail for use in the context of the multimedia directory plugin and other plugins where appropriate. It operates on files with the *.sldprt, *.sldasm, and *.slddrw extensions.

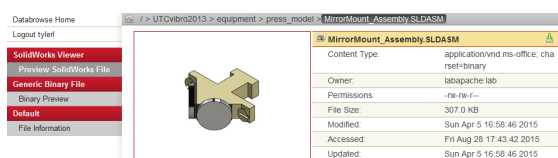

7.23 Specimen Management Plugin

The Specimen Management plugin provides an interface for storing and tracking information about specimens. Operating on XML files containing the *.sdb extension, it can produce a visual display of the details of a specimen, including identification, geometry, provenance information, images, bar code labels, etc. Additionally, utilizing the Axel XML JavaScript library (<http://ssire.github.io/axel/>), this plugin can also provide a web-based interface to edit the details of a specimen on an interactive form. The structure of the sdb files are outside of the context of this document. Templates that control the display and editing of specimen files can be modified, enabling the introduction of new parameters. This plugin is also capable of combining multiple sources of data about a specimen into one unified representation, highly useful both on the web and in the context of the Databrowse library. This capability is presently being used to enable specimens to be placed into specimen groups – sharing all of the parameters from the group and thus limiting the unnecessary repetition of data across many files.

The screenshot displays the 'Specimen Details' form for specimen C00-ETC2B-092L. The interface includes a sidebar with navigation options like 'Specimen', 'Edit Specimen Data', and 'View Specimen Data'. The main content area is divided into several sections: 'Specimen Details' (ID, Group, Support Files), 'Geometric Properties' (Component, Shape, Dimension Type, Coordinate Sys, Dimensions), 'Physical Properties' (Component, Material), 'Flaw Parameters' (Index, Flaw Type, Flaw Shape, Flaw Size, Coordinate, Frame, Location), 'Provenance' (Manufactured By, Source, In Charge), 'Specimen Markings' (Identifier Marks), 'Specimen Orientation' (Face Definitions), and 'Supporting Files' (a list of image files). A 'Specimen Barcode' section on the right shows a QR code and a 'Print Barcode on Label Printer' button.

7.24 Specimen Directory Plugin

The Specimen Directory plugin is a complementary tool to the Specimen Management plugin, overriding the directory plugin on directories that contain *.sdb files. This plugin displays a tabular view of specimen parameters, enabling the quick searching, filtering, and sorting of entire sets of specimen data. This plugin also provides easy access needed to create new *.sdb files in the current location.

The screenshot shows the 'Specimen Database' table with columns for ID, Groups, Material, Shape, Dimension 1, Dimension 2, Dimension 3, Crack Size, and Crack Location. The table lists several specimens, including C00-ETC2B-109C, C00-ETC2B-109B, C00-ETC2B-109A, C00-ETC2B-109D, C00-ETC2B-109E, C00-ETC2B-109F, C00-ETC2B-109G, C00-ETC2B-109H, C00-ETC2B-109I, C00-ETC2B-109J, C00-ETC2B-109K, C00-ETC2B-109L, C00-ETC2B-109M, C00-ETC2B-109N, C00-ETC2B-109O, C00-ETC2B-109P, C00-ETC2B-109Q, C00-ETC2B-109R, C00-ETC2B-109S, C00-ETC2B-109T, C00-ETC2B-109U, C00-ETC2B-109V, C00-ETC2B-109W, C00-ETC2B-109X, C00-ETC2B-109Y, C00-ETC2B-109Z, C00-ETC2B-109AA, C00-ETC2B-109AB, C00-ETC2B-109AC, C00-ETC2B-109AD, C00-ETC2B-109AE, C00-ETC2B-109AF, C00-ETC2B-109AG, C00-ETC2B-109AH, C00-ETC2B-109AI, C00-ETC2B-109AJ, C00-ETC2B-109AK, C00-ETC2B-109AL, C00-ETC2B-109AM, C00-ETC2B-109AN, C00-ETC2B-109AO, C00-ETC2B-109AP, C00-ETC2B-109AQ, C00-ETC2B-109AR, C00-ETC2B-109AS, C00-ETC2B-109AT, C00-ETC2B-109AU, C00-ETC2B-109AV, C00-ETC2B-109AW, C00-ETC2B-109AX, C00-ETC2B-109AY, C00-ETC2B-109AZ, C00-ETC2B-109BA, C00-ETC2B-109BB, C00-ETC2B-109BC, C00-ETC2B-109BD, C00-ETC2B-109BE, C00-ETC2B-109BF, C00-ETC2B-109BG, C00-ETC2B-109BH, C00-ETC2B-109BI, C00-ETC2B-109BJ, C00-ETC2B-109BK, C00-ETC2B-109BL, C00-ETC2B-109BM, C00-ETC2B-109BN, C00-ETC2B-109BO, C00-ETC2B-109BP, C00-ETC2B-109BQ, C00-ETC2B-109BR, C00-ETC2B-109BS, C00-ETC2B-109BT, C00-ETC2B-109BU, C00-ETC2B-109BV, C00-ETC2B-109BW, C00-ETC2B-109BX, C00-ETC2B-109BY, C00-ETC2B-109BZ, C00-ETC2B-109CA, C00-ETC2B-109CB, C00-ETC2B-109CC, C00-ETC2B-109CD, C00-ETC2B-109CE, C00-ETC2B-109CF, C00-ETC2B-109CG, C00-ETC2B-109CH, C00-ETC2B-109CI, C00-ETC2B-109CJ, C00-ETC2B-109CK, C00-ETC2B-109CL, C00-ETC2B-109CM, C00-ETC2B-109CN, C00-ETC2B-109CO, C00-ETC2B-109CP, C00-ETC2B-109CQ, C00-ETC2B-109CR, C00-ETC2B-109CS, C00-ETC2B-109CT, C00-ETC2B-109CU, C00-ETC2B-109CV, C00-ETC2B-109CW, C00-ETC2B-109CX, C00-ETC2B-109CY, C00-ETC2B-109CZ, C00-ETC2B-109DA, C00-ETC2B-109DB, C00-ETC2B-109DC, C00-ETC2B-109DD, C00-ETC2B-109DE, C00-ETC2B-109DF, C00-ETC2B-109DG, C00-ETC2B-109DH, C00-ETC2B-109DI, C00-ETC2B-109DJ, C00-ETC2B-109DK, C00-ETC2B-109DL, C00-ETC2B-109DM, C00-ETC2B-109DN, C00-ETC2B-109DO, C00-ETC2B-109DP, C00-ETC2B-109DQ, C00-ETC2B-109DR, C00-ETC2B-109DS, C00-ETC2B-109DT, C00-ETC2B-109DU, C00-ETC2B-109DV, C00-ETC2B-109DW, C00-ETC2B-109DX, C00-ETC2B-109DY, C00-ETC2B-109DZ, C00-ETC2B-109EA, C00-ETC2B-109EB, C00-ETC2B-109EC, C00-ETC2B-109ED, C00-ETC2B-109EE, C00-ETC2B-109EF, C00-ETC2B-109EG, C00-ETC2B-109EH, C00-ETC2B-109EI, C00-ETC2B-109EJ, C00-ETC2B-109EK, C00-ETC2B-109EL, C00-ETC2B-109EM, C00-ETC2B-109EN, C00-ETC2B-109EO, C00-ETC2B-109EP, C00-ETC2B-109EQ, C00-ETC2B-109ER, C00-ETC2B-109ES, C00-ETC2B-109ET, C00-ETC2B-109EU, C00-ETC2B-109EV, C00-ETC2B-109EW, C00-ETC2B-109EX, C00-ETC2B-109EY, C00-ETC2B-109EZ, C00-ETC2B-109FA, C00-ETC2B-109FB, C00-ETC2B-109FC, C00-ETC2B-109FD, C00-ETC2B-109FE, C00-ETC2B-109FF, C00-ETC2B-109FG, C00-ETC2B-109FH, C00-ETC2B-109FI, C00-ETC2B-109FJ, C00-ETC2B-109FK, C00-ETC2B-109FL, C00-ETC2B-109FM, C00-ETC2B-109FN, C00-ETC2B-109FO, C00-ETC2B-109FP, C00-ETC2B-109FQ, C00-ETC2B-109FR, C00-ETC2B-109FS, C00-ETC2B-109FT, C00-ETC2B-109FU, C00-ETC2B-109FV, C00-ETC2B-109FW, C00-ETC2B-109FX, C00-ETC2B-109FY, C00-ETC2B-109FZ, C00-ETC2B-109GA, C00-ETC2B-109GB, C00-ETC2B-109GC, C00-ETC2B-109GD, C00-ETC2B-109GE, C00-ETC2B-109GF, C00-ETC2B-109GG, C00-ETC2B-109GH, C00-ETC2B-109GI, C00-ETC2B-109GJ, C00-ETC2B-109GK, C00-ETC2B-109GL, C00-ETC2B-109GM, C00-ETC2B-109GN, C00-ETC2B-109GO, C00-ETC2B-109GP, C00-ETC2B-109GQ, C00-ETC2B-109GR, C00-ETC2B-109GS, C00-ETC2B-109GT, C00-ETC2B-109GU, C00-ETC2B-109GV, C00-ETC2B-109GW, C00-ETC2B-109GX, C00-ETC2B-109GY, C00-ETC2B-109GZ, C00-ETC2B-109HA, C00-ETC2B-109HB, C00-ETC2B-109HC, C00-ETC2B-109HD, C00-ETC2B-109HE, C00-ETC2B-109HF, C00-ETC2B-109HG, C00-ETC2B-109HH, C00-ETC2B-109HI, C00-ETC2B-109HJ, C00-ETC2B-109HK, C00-ETC2B-109HL, C00-ETC2B-109HM, C00-ETC2B-109HN, C00-ETC2B-109HO, C00-ETC2B-109HP, C00-ETC2B-109HQ, C00-ETC2B-109HR, C00-ETC2B-109HS, C00-ETC2B-109HT, C00-ETC2B-109HU, C00-ETC2B-109HV, C00-ETC2B-109HW, C00-ETC2B-109HX, C00-ETC2B-109HY, C00-ETC2B-109HZ, C00-ETC2B-109IA, C00-ETC2B-109IB, C00-ETC2B-109IC, C00-ETC2B-109ID, C00-ETC2B-109IE, C00-ETC2B-109IF, C00-ETC2B-109IG, C00-ETC2B-109IH, C00-ETC2B-109II, C00-ETC2B-109IJ, C00-ETC2B-109IK, C00-ETC2B-109IL, C00-ETC2B-109IM, C00-ETC2B-109IN, C00-ETC2B-109IO, C00-ETC2B-109IP, C00-ETC2B-109IQ, C00-ETC2B-109IR, C00-ETC2B-109IS, C00-ETC2B-109IT, C00-ETC2B-109IU, C00-ETC2B-109IV, C00-ETC2B-109IW, C00-ETC2B-109IX, C00-ETC2B-109IY, C00-ETC2B-109IZ, C00-ETC2B-109JA, C00-ETC2B-109JB, C00-ETC2B-109JC, C00-ETC2B-109JD, C00-ETC2B-109JE, C00-ETC2B-109JF, C00-ETC2B-109JG, C00-ETC2B-109JH, C00-ETC2B-109JI, C00-ETC2B-109JJ, C00-ETC2B-109JK, C00-ETC2B-109JL, C00-ETC2B-109JM, C00-ETC2B-109JN, C00-ETC2B-109JO, C00-ETC2B-109JP, C00-ETC2B-109JQ, C00-ETC2B-109JR, C00-ETC2B-109JS, C00-ETC2B-109JT, C00-ETC2B-109JU, C00-ETC2B-109JV, C00-ETC2B-109JW, C00-ETC2B-109JX, C00-ETC2B-109JY, C00-ETC2B-109JZ, C00-ETC2B-109KA, C00-ETC2B-109KB, C00-ETC2B-109KC, C00-ETC2B-109KD, C00-ETC2B-109KE, C00-ETC2B-109KF, C00-ETC2B-109KG, C00-ETC2B-109KH, C00-ETC2B-109KI, C00-ETC2B-109KJ, C00-ETC2B-109KK, C00-ETC2B-109KL, C00-ETC2B-109KM, C00-ETC2B-109KN, C00-ETC2B-109KO, C00-ETC2B-109KP, C00-ETC2B-109KQ, C00-ETC2B-109KR, C00-ETC2B-109KS, C00-ETC2B-109KT, C00-ETC2B-109KU, C00-ETC2B-109KV, C00-ETC2B-109KW, C00-ETC2B-109KX, C00-ETC2B-109KY, C00-ETC2B-109KZ, C00-ETC2B-109LA, C00-ETC2B-109LB, C00-ETC2B-109LC, C00-ETC2B-109LD, C00-ETC2B-109LE, C00-ETC2B-109LF, C00-ETC2B-109LG, C00-ETC2B-109LH, C00-ETC2B-109LI, C00-ETC2B-109LJ, C00-ETC2B-109LK, C00-ETC2B-109LL, C00-ETC2B-109LM, C00-ETC2B-109LN, C00-ETC2B-109LO, C00-ETC2B-109LP, C00-ETC2B-109LQ, C00-ETC2B-109LR, C00-ETC2B-109LS, C00-ETC2B-109LT, C00-ETC2B-109LU, C00-ETC2B-109LV, C00-ETC2B-109LW, C00-ETC2B-109LX, C00-ETC2B-109LY, C00-ETC2B-109LZ, C00-ETC2B-109MA, C00-ETC2B-109MB, C00-ETC2B-109MC, C00-ETC2B-109MD, C00-ETC2B-109ME, C00-ETC2B-109MF, C00-ETC2B-109MG, C00-ETC2B-109MH, C00-ETC2B-109MI, C00-ETC2B-109MJ, C00-ETC2B-109MK, C00-ETC2B-109ML, C00-ETC2B-109MN, C00-ETC2B-109MO, C00-ETC2B-109MP, C00-ETC2B-109MQ, C00-ETC2B-109MR, C00-ETC2B-109MS, C00-ETC2B-109MT, C00-ETC2B-109MU, C00-ETC2B-109MV, C00-ETC2B-109MW, C00-ETC2B-109MX, C00-ETC2B-109MY, C00-ETC2B-109MZ, C00-ETC2B-109NA, C00-ETC2B-109NB, C00-ETC2B-109NC, C00-ETC2B-109ND, C00-ETC2B-109NE, C00-ETC2B-109NF, C00-ETC2B-109NG, C00-ETC2B-109NH, C00-ETC2B-109NI, C00-ETC2B-109NJ, C00-ETC2B-109NK, C00-ETC2B-109NL, C00-ETC2B-109NM, C00-ETC2B-109NN, C00-ETC2B-109NO, C00-ETC2B-109NP, C00-ETC2B-109NQ, C00-ETC2B-109NR, C00-ETC2B-109NS, C00-ETC2B-109NT, C00-ETC2B-109NU, C00-ETC2B-109NV, C00-ETC2B-109NW, C00-ETC2B-109NX, C00-ETC2B-109NY, C00-ETC2B-109NZ, C00-ETC2B-109OA, C00-ETC2B-109OB, C00-ETC2B-109OC, C00-ETC2B-109OD, C00-ETC2B-109OE, C00-ETC2B-109OF, C00-ETC2B-109OG, C00-ETC2B-109OH, C00-ETC2B-109OI, C00-ETC2B-109OJ, C00-ETC2B-109OK, C00-ETC2B-109OL, C00-ETC2B-109OM, C00-ETC2B-109ON, C00-ETC2B-109OO, C00-ETC2B-109OP, C00-ETC2B-109OQ, C00-ETC2B-109OR, C00-ETC2B-109OS, C00-ETC2B-109OT, C00-ETC2B-109OU, C00-ETC2B-109OV, C00-ETC2B-109OW, C00-ETC2B-109OX, C00-ETC2B-109OY, C00-ETC2B-109OZ, C00-ETC2B-109PA, C00-ETC2B-109PB, C00-ETC2B-109PC, C00-ETC2B-109PD, C00-ETC2B-109PE, C00-ETC2B-109PF, C00-ETC2B-109PG, C00-ETC2B-109PH, C00-ETC2B-109PI, C00-ETC2B-109PJ, C00-ETC2B-109PK, C00-ETC2B-109PL, C00-ETC2B-109PM, C00-ETC2B-109PN, C00-ETC2B-109PO, C00-ETC2B-109PP, C00-ETC2B-109PQ, C00-ETC2B-109PR, C00-ETC2B-109PS, C00-ETC2B-109PT, C00-ETC2B-109PU, C00-ETC2B-109PV, C00-ETC2B-109PW, C00-ETC2B-109PX, C00-ETC2B-109PY, C00-ETC2B-109PZ, C00-ETC2B-109QA, C00-ETC2B-109QB, C00-ETC2B-109QC, C00-ETC2B-109QD, C00-ETC2B-109QE, C00-ETC2B-109QF, C00-ETC2B-109QG, C00-ETC2B-109QH, C00-ETC2B-109QI, C00-ETC2B-109QJ, C00-ETC2B-109QK, C00-ETC2B-109QL, C00-ETC2B-109QM, C00-ETC2B-109QN, C00-ETC2B-109QO, C00-ETC2B-109QP, C00-ETC2B-109QQ, C00-ETC2B-109QR, C00-ETC2B-109QS, C00-ETC2B-109QT, C00-ETC2B-109QU, C00-ETC2B-109QV, C00-ETC2B-109QW, C00-ETC2B-109QX, C00-ETC2B-109QY, C00-ETC2B-109QZ, C00-ETC2B-109RA, C00-ETC2B-109RB, C00-ETC2B-109RC, C00-ETC2B-109RD, C00-ETC2B-109RE, C00-ETC2B-109RF, C00-ETC2B-109RG, C00-ETC2B-109RH, C00-ETC2B-109RI, C00-ETC2B-109RJ, C00-ETC2B-109RK, C00-ETC2B-109RL, C00-ETC2B-109RM, C00-ETC2B-109RN, C00-ETC2B-109RO, C00-ETC2B-109RP, C00-ETC2B-109RQ, C00-ETC2B-109RR, C00-ETC2B-109RS, C00-ETC2B-109RT, C00-ETC2B-109RU, C00-ETC2B-109RV, C00-ETC2B-109RW, C00-ETC2B-109RX, C00-ETC2B-109RY, C00-ETC2B-109RZ, C00-ETC2B-109SA, C00-ETC2B-109SB, C00-ETC2B-109SC, C00-ETC2B-109SD, C00-ETC2B-109SE, C00-ETC2B-109SF, C00-ETC2B-109SG, C00-ETC2B-109SH, C00-ETC2B-109SI, C00-ETC2B-109SJ, C00-ETC2B-109SK, C00-ETC2B-109SL, C00-ETC2B-109SM, C00-ETC2B-109SN, C00-ETC2B-109SO, C00-ETC2B-109SP, C00-ETC2B-109SQ, C00-ETC2B-109SR, C00-ETC2B-109SS, C00-ETC2B-109ST, C00-ETC2B-109SU, C00-ETC2B-109SV, C00-ETC2B-109SW, C00-ETC2B-109SX, C00-ETC2B-109SY, C00-ETC2B-109SZ, C00-ETC2B-109TA, C00-ETC2B-109TB, C00-ETC2B-109TC, C00-ETC2B-109TD, C00-ETC2B-109TE, C00-ETC2B-109TF, C00-ETC2B-109TG, C00-ETC2B-109TH, C00-ETC2B-109TI, C00-ETC2B-109TJ, C00-ETC2B-109TK, C00-ETC2B-109TL, C00-ETC2B-109TM, C00-ETC2B-109TN, C00-ETC2B-109TO, C00-ETC2B-109TP, C00-ETC2B-109TQ, C00-ETC2B-109TR, C00-ETC2B-109TS, C00-ETC2B-109TT, C00-ETC2B-109TU, C00-ETC2B-109TV, C00-ETC2B-109TW, C00-ETC2B-109TX, C00-ETC2B-109TY, C00-ETC2B-109TZ, C00-ETC2B-109UA, C00-ETC2B-109UB, C00-ETC2B-109UC, C00-ETC2B-109UD, C00-ETC2B-109UE, C00-ETC2B-109UF, C00-ETC2B-109UG, C00-ETC2B-109UH, C00-ETC2B-109UI, C00-ETC2B-109UJ, C00-ETC2B-109UK, C00-ETC2B-109UL, C00-ETC2B-109UM, C00-ETC2B-109UN, C00-ETC2B-109UO, C00-ETC2B-109UP, C00-ETC2B-109UQ, C00-ETC2B-109UR, C00-ETC2B-109US, C00-ETC2B-109UT, C00-ETC2B-109UU, C00-ETC2B-109UV, C00-ETC2B-109UW, C00-ETC2B-109UX, C00-ETC2B-109UY, C00-ETC2B-109UZ, C00-ETC2B-109VA, C00-ETC2B-109VB, C00-ETC2B-109VC, C00-ETC2B-109VD, C00-ETC2B-109VE, C00-ETC2B-109VF, C00-ETC2B-109VG, C00-ETC2B-109VH, C00-ETC2B-109VI, C00-ETC2B-109VJ, C00-ETC2B-109VK, C00-ETC2B-109VL, C00-ETC2B-109VM, C00-ETC2B-109VN, C00-ETC2B-109VO, C00-ETC2B-109VP, C00-ETC2B-109VQ, C00-ETC2B-109VR, C00-ETC2B-109VS, C00-ETC2B-109VT, C00-ETC2B-109VU, C00-ETC2B-109VV, C00-ETC2B-109VW, C00-ETC2B-109VX, C00-ETC2B-109VY, C00-ETC2B-109VZ, C00-ETC2B-109WA, C00-ETC2B-109WB, C00-ETC2B-109WC, C00-ETC2B-109WD, C00-ETC2B-109WE, C00-ETC2B-109WF, C00-ETC2B-109WG, C00-ETC2B-109WH, C00-ETC2B-109WI, C00-ETC2B-109WJ, C00-ETC2B-109WK, C00-ETC2B-109WL, C00-ETC2B-109WM, C00-ETC2B-109WN, C00-ETC2B-109WO, C00-ETC2B-109WP, C00-ETC2B-109WQ, C00-ETC2B-109WR, C00-ETC2B-109WS, C00-ETC2B-109WT, C00-ETC2B-109WU, C00-ETC2B-109WV, C00-ETC2B-109WW, C00-ETC2B-109WX, C00-ETC2B-109WY, C00-ETC2B-109WZ, C00-ETC2B-109XA, C00-ETC2B-109XB, C00-ETC2B-109XC, C00-ETC2B-109XD, C00-ETC2B-109XE, C00-ETC2B-109XF, C00-ETC2B-109XG, C00-ETC2B-109XH, C00-ETC2B-109XI, C00-ETC2B-109XJ, C00-ETC2B-109XK, C00-ETC2B-109XL, C00-ETC2B-109XM, C00-ETC2B-109XN, C00-ETC2B-109XO, C00-ETC2B-109XP, C00-ETC2B-109XQ, C00-ETC2B-109XR, C00-ETC2B-109XS, C00-ETC2B-109XT, C00-ETC2B-109XU, C00-ETC2B-109XV, C00-ETC2B-109XW, C00-ETC2B-109XX, C00-ETC2B-109XY, C00-ETC2B-109XZ, C00-ETC2B-109YA, C00-ETC2B-109YB, C00-ETC2B-109YC, C00-ETC2B-109YD, C00-ETC2B-109YE, C00-ETC2B-109YF, C00-ETC2B-109YG, C00-ETC2B-109YH, C00-ETC2B-109YI, C00-ETC2B-109YJ, C00-ETC2B-109YK, C00-ETC2B-109YL, C00-ETC2B-109YM, C00-ETC2B-109YN, C00-ETC2B-109YO, C00-ETC2B-109YP, C00-ETC2B-109YQ, C00-ETC2B-109YR, C00-ETC2B-109YS, C00-ETC2B-109YT, C00-ETC2B-109YU, C00-ETC2B-109YV, C00-ETC2B-109YW, C00-ETC2B-109YX, C00-ETC2B-109YY, C00-ETC2B-109YZ, C00-ETC2B-109ZA, C00-ETC2B-109ZB, C00-ETC2B-109ZC, C00-ETC2B-109ZD, C00-ETC2B-109ZE, C00-ETC2B-109ZF, C00-ETC2B-109ZG, C00-ETC2B-109ZH, C00-ETC2B-109ZI, C00-ETC2B-109ZJ, C00-ETC2B-109ZK, C00-ETC2B-109ZL, C00-ETC2B-109ZM, C00-ETC2B-109ZN, C00-ETC2B-109ZO, C00-ETC2B-109ZP, C00-ETC2B-109ZQ, C00-ETC2B-109ZR, C00-ETC2B-109ZS, C00-ETC2B-109ZT, C00-ETC2B-109ZU, C00-ETC2B-109ZV, C00-ETC2B-109ZW, C00-ETC2B-109ZX, C00-ETC2B-109ZY, C00-ETC2B-109ZZ.

7.25 Specimen Group Management Plugin

The Specimen Group Management plugin is a complementary tool to the Specimen Management plugin, providing similar functionality, but acting in the context of *.sdg files. The schema defining the *.sdg file is outside of the context of this document; however, it is similar in structure to the *.sdb file.



The screenshot displays the 'Specimen Group Details' form for group UTCA. The interface includes a sidebar with navigation options like 'Specimen Group', 'Edit Specimen Group', and 'View Specimen Group'. The main content area is divided into several sections: 'Specimen Group Details' (Group ID, Description, Tags, Support Files), 'Physical Properties' (Component, Material, Young's Modulus, Poisson's Ratio, Shear Modulus, Density, Yield Strength), 'Provenance' (Manufactured By, Source, In Charge), 'Specimen Markings' (Identifier Marks), 'Specimen Orientation' (Face Definitions), and 'Supporting Files' (a list of PDF files). The 'Physical Properties' section includes detailed information about the material (Ti-6-4) and its mechanical properties, such as Young's Modulus (121.45 GPa) and Poisson's Ratio (0.342).

7.26 SVG Viewer

The SVG viewer plugin adds the necessary support to ensure that *.svg vector graphics image files can be displayed on all web browsers, producing an image thumbnail and preview if necessary. It also produces thumbnails for the multimedia directory plugin and in other contexts as needed.


7.27 Transducer Management Plugin

The Transducer Management plugin is similar in functionality to the Specimen Management plugin, but designed in the context of managing parameters and other information associated with ultrasonic transducers. The plugin operates XML files with the extension *.tdb. The structure of these files is outside of the context of this document.

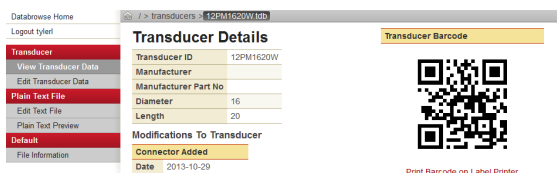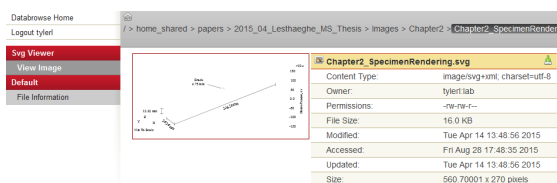


7.28 Transducer Directory Plugin

The Transducer Directory plugin is similar in functionality to the Specimen Directory plugin, overriding the directory plugin on directories that contain *.tdb files, and providing a tabular view of the data contained within, enabling rapid searching, sorting, and filtering.

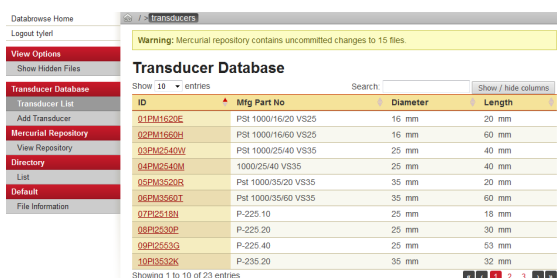


7.29 Trigger Log Plugin

The Trigger Log plugin is similar in functionality to the Transducer Management plugin, but designed in the context of tracking the usage of ultrasonic transducers. An interface with data acquisition software ensures that an XML file with the *.tgl extension is recorded to ensure accurate tracking of all experimental triggers. This plugin will display the data from such a file in a tabular form, enabling rapid searching, sorting, and filtering. The structure of the *.tgl file is outside of the context of this document.

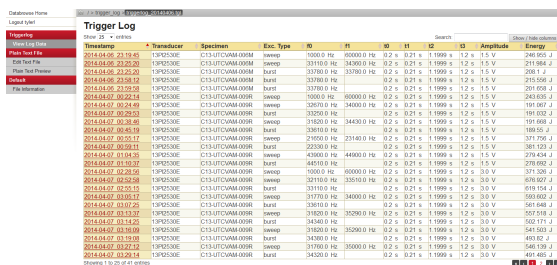

7.30 Trigger Log Directory Plugin

The Trigger Log Directory plugin functions similarly to the Trigger Log plugin, but enabling an entire directory of *.tgl files to be combined into one single representation. This representation is displayed in tabular form, enabling rapid searching, sorting and filtering. This plugin also provides the capability of exporting this data to a CSV file for further analysis.

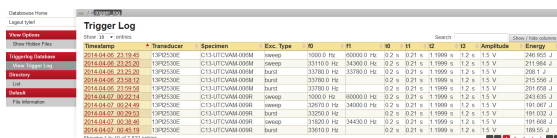

7.31 Web Page Viewer

The web page viewer plugin will display the contents of an HTML file within the Databrowse interface. Combined with the plain text editor plugin, this plugin can provide a convenient mechanism for documentation while being able to pull together all of the other resources available inside the context of Databrowse.

Databrowse Home

Legend: 1x1pt

Webpage Viewer

View Webpage

Plain Text File

Edit Text File

Plain Text Preview

Default

File Information

go / > UTCvibro2013 > specimen_geometry_selection > analysis > SpecimenGeometryRationale.html

Specimen Geometry Choice Rationale

This document outlines the details behind the choice of geometry used in the UTCIAFRL Vibrothermography Modeling Project. Several factors played into the choice of geometry used for experimental specimens, which are outlined below. The chosen geometry are as follows:

- Inconel 718: 275 mm x 25.4 mm x 12.7 mm
- Ti 6-4: 245 mm x 25.4 mm x 12.7 mm

Sample Thickness Selection

We ultimately desire to minimize the required thickness, as additional thickness results in higher load required to apply the desired bending stresses (approx. 50-250 MPa) and adds additional material cost. For the lower bound, we need to consider that we desire to grow long cracks (approx. 7-12 mm) in the specimen. Samples that are 1/2 inch (12.7 mm) thick provide the necessary thickness to grow cracks of this size.

Sample Width Selection

The width of the samples were chosen after doing several rounds of motion scanning on the specimens. During this process, the sample was excited using a frequency sweep between 100 Hz and 50 kHz for 1 second. This was repeated several hundred times (depending on the geometry of the test specimen) with a laser vibrometer pointed at various points on the sample. This allowed us to identify the resonance modes present in a particular geometry. Three sample widths were tested:

- 1 inch (10 inches long)
- 1.5 inches (10 inches long)
- 2 inches at (6 inches long)

The ideal specimen width will allow us to easily excite pure flexural bending modes. As sample width increases, the ability to do this diminishes as more torsional and hybrid modes appear, and appear closer in frequency to the pure flexural modes we desire. It is also worth noting that 1 inch serves as a good lower bound, considering the length of cracks we desire to grow.

Our vibration validation experiments show that the 1 inch (25.4 mm) wide samples generally provide pure flexure in the 7th, 9th, and 11th modes while minimizing the appearance of torsional and other hybrid effects near the frequencies required to excite these modes. These modes are ideal for our experiments as they produce high strains at the locations where cracks will be produced. The 1.5 inch sample was considered as a possible choice, since they would provide the benefit of being able to grow longer cracks. However, their performance did not meet expectations, as it was somewhat difficult to identify pure flexural modes in these specimens. Identifying pure flexural modes in the 2 inch wide specimens was extremely difficult.

It is also important to note that 1 inch wide samples generate higher heating index values (strain²/frequency), which is ideal for our experiments as well.

Both of these points are illustrated in the plots below.

Heat Index vs Specimen Width (Inconel 718)

35

8 License and Third-Party Components

Databrowse is distributed under the terms of the BSD-3 clause license.

Databrowse is packaged with the following third party components:

Portions of this software are adapted from gnome-icon-theme-3.7.4:

Copyright © 2005-2013 The GNOME Project

License: LGPL v3 or CC-BY-SA v3.0

Portions of this software are adapted from jQuery-1.9.1:

Copyright © 2013 jQuery Foundation and other contributors

License: MIT

Portions of this software are adapted from DataTables-1.9.4:

Copyright © 2008-2010 Allan Jardine

License: GPL v2 or BSD-3

Portions of this software are adapted from hurry.filesize-0.9:

Copyright © 2009 Martijn Faassen, Startifact

License: ZPL v2.1

Portions of this software are adapted from prettify-4-Mar-2013:

Copyright © 2011 Mike Samuel et al

License: Apache v2.0

Portions of this software are adapted from jQuery-File-Upload-7.4.1-jquery.ui:

Copyright © 2013 Sebastian Tschan

License: MIT

Portions of this software are adapted from axel-1.1.2-beta:

Copyright © 2008 - 2010 Staphane Sire

License: LGPL v2.1 or newer

Portions of this software are adapted from images2gif.py:

Copyright © 2012 Almar Klein, Ant1, Marius van Voorden

License: BSD-3

Portions of this software are adapted from exif-py-1.2.0:

Copyright © 2002-2007 Gene Cash

Copyright © 2007-2013 Ianar Svi and contributors

License: BSD-3

8.1 BSD 3-Clause License

Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

1. Redistributions of source code must retain the above copyright notice, this list of conditions and the following disclaimer.
2. Redistributions in binary form must reproduce the above copyright notice, this list of conditions and the following disclaimer in the documentation and/or other materials provided with the distribution.
3. The name of the author may not be used to endorse or promote products derived from this software without specific prior written permission.

THIS SOFTWARE IS PROVIDED BY THE AUTHOR “AS IS” AND ANY EXPRESS OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE AUTHOR BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

8.2 GNU General Public License v3

GNU GENERAL PUBLIC LICENSE

Version 3, 29 June 2007

Copyright © 2007 Free Software Foundation, Inc. <http://fsf.org/>

Everyone is permitted to copy and distribute verbatim copies of this license document, but changing it is not allowed.

Preamble

The GNU General Public License is a free, copyleft license for software and other kinds of works.

The licenses for most software and other practical works are designed to take away your freedom to share and change the works. By contrast, the GNU General Public License is intended to guarantee your freedom to share and change all versions of a program—to make sure it remains free software for all its users. We, the Free Software Foundation, use the GNU General Public License for most of our software; it applies also to any other work released this way by its authors. You can apply it to your programs, too.

When we speak of free software, we are referring to freedom, not price. Our General Public Licenses are designed to make sure that you have the freedom to distribute copies of free software (and charge for them if you wish), that you receive source code or can get it if you want it, that you can change the software or use pieces of it in new free programs, and that you know you can do these things.

To protect your rights, we need to prevent others from denying you these rights or asking you to surrender the rights. Therefore, you have certain responsibilities if you distribute copies of the software, or if you modify it: responsibilities to respect the freedom of others.

For example, if you distribute copies of such a program, whether gratis or for a fee, you must pass on to the recipients the same freedoms that you received. You must make sure that they, too, receive or can get the source code. And you must show them these terms so they know their rights.

Developers that use the GNU GPL protect your rights with two steps: (1) assert copyright on the software, and (2) offer you this License giving you legal permission to copy, distribute and/or modify it.

For the developers' and authors' protection, the GPL clearly explains that there is no warranty for this free software. For both users' and authors' sake, the GPL requires that modified versions be marked as changed, so that their problems will not be attributed erroneously to authors of previous versions.

Some devices are designed to deny users access to install or run modified versions of the software inside them, although the manufacturer can do so. This is fundamentally incompatible with the aim of protecting users' freedom to change the software. The systematic pattern of such abuse occurs in the area of products for individuals to use, which is precisely where it is most unacceptable. Therefore, we have designed this version of the GPL to prohibit the practice for those products. If such problems arise substantially in other domains, we stand ready to extend this provision to those domains in future versions of the GPL, as needed to protect the freedom of users.

Finally, every program is threatened constantly by software patents. States should not allow patents to restrict development and use of software on general-purpose computers, but in those that do, we wish to avoid the special danger that patents applied to a free program could make it effectively proprietary. To prevent this, the GPL assures that patents cannot be used to render the program non-free.

The precise terms and conditions for copying, distribution and modification follow.

TERMS AND CONDITIONS

0. Definitions.

“This License” refers to version 3 of the GNU General Public License.

“Copyright” also means copyright-like laws that apply to other kinds of works, such as semiconductor masks.

“The Program” refers to any copyrightable work licensed under this License. Each licensee is addressed as “you”. “Licensees” and “recipients” may be individuals or organizations.

To “modify” a work means to copy from or adapt all or part of the work in a fashion requiring copyright permission, other than the making of an exact copy. The resulting work is called a “modified version” of the earlier work or a work “based on” the earlier work.

A “covered work” means either the unmodified Program or a work based on the Program.

To “propagate” a work means to do anything with it that, without permission, would make you directly or secondarily liable for infringement under applicable copyright law, except executing it on a computer or modifying a private copy. Propagation includes copying, distribution (with or without modification), making available to the public, and in some countries other activities as well.

To “convey” a work means any kind of propagation that enables other parties to make or receive copies. Mere interaction with a user through a computer network, with no transfer of a copy, is not conveying.

An interactive user interface displays “Appropriate Legal Notices” to the extent that it includes a convenient and prominently visible feature that (1) displays an appropriate copyright notice, and (2) tells the user that there is no warranty for the work (except to the extent that warranties are provided), that licensees may convey the work under this License, and how to view a copy of this License. If the interface presents a list of user commands or options, such as a menu, a prominent item in the list meets this criterion.

1. Source Code.

The “source code” for a work means the preferred form of the work for making modifications to it. “Object code” means any non-source form of a work.

A “Standard Interface” means an interface that either is an official standard defined by a recognized standards body, or, in the case of interfaces specified for a particular programming language, one that is widely used among developers working in that language.

The “System Libraries” of an executable work include anything, other than the work as a whole, that (a) is included in the normal form of packaging a Major Component, but which is not part of that Major Component, and (b) serves only to enable use of the work with that Major Component, or to implement a Standard Interface for which an implementation is available to the public in source code form. A “Major Component”, in this context, means a major essential component (kernel, window system, and so on) of the specific operating system (if any) on which the executable work runs, or a compiler used to produce the work, or an object code interpreter used to run it.

The “Corresponding Source” for a work in object code form means all the source code needed to generate, install, and (for an executable work) run the object code and to modify the work, including scripts to control those activities. However, it does not include the work’s System Libraries, or general-purpose tools or generally available free programs which are used unmodified in performing those activities but which are not part of the work. For example, Corresponding Source includes interface definition files associated with source files for the work, and the source code for shared libraries and dynamically linked subprograms that the work is specifically designed to require, such as by intimate data communication or control flow between those subprograms and other parts of the work.

The Corresponding Source need not include anything that users can regenerate automatically from other parts of the Corresponding Source.

The Corresponding Source for a work in source code form is that same work.

2. Basic Permissions.

All rights granted under this License are granted for the term of copyright on the Program, and are irrevocable provided the stated conditions are met. This License explicitly affirms your unlimited permission to run the unmodified Program. The output from running a covered work is covered by this License only if the output, given its content, constitutes a covered work. This License acknowledges your rights of fair use or other equivalent, as provided by copyright law.

You may make, run and propagate covered works that you do not convey, without conditions so long as your license otherwise remains in force. You may convey covered works to others for the sole purpose of having them make modifications exclusively for you, or provide you with facilities for running those works, provided that you comply with the terms of this License in conveying all material for which you do not control copyright. Those thus making or running the covered works for you must do so exclusively on your behalf, under your direction and control, on terms that prohibit them from making any copies of your copyrighted material outside their relationship with you.

Conveying under any other circumstances is permitted solely under the conditions stated below. Sub-licensing is not allowed; section 10 makes it unnecessary.

3. Protecting Users' Legal Rights From Anti-Circumvention Law.

No covered work shall be deemed part of an effective technological measure under any applicable law fulfilling obligations under article 11 of the WIPO copyright treaty adopted on 20 December 1996, or similar laws prohibiting or restricting circumvention of such measures.

When you convey a covered work, you waive any legal power to forbid circumvention of technological measures to the extent such circumvention is effected by exercising rights under this License with respect to the covered work, and you disclaim any intention to limit operation or modification of the work as a means of enforcing, against the work's users, your or third parties' legal rights to forbid circumvention of technological measures.

4. Conveying Verbatim Copies.

You may convey verbatim copies of the Program's source code as you receive it, in any medium, provided that you conspicuously and appropriately publish on each copy an appropriate copyright notice; keep intact all notices stating that this License and any non-permissive terms added in accord with section 7 apply to the code; keep intact all notices of the absence of any warranty; and give all recipients a copy of this License along with the Program.

You may charge any price or no price for each copy that you convey, and you may offer support or warranty protection for a fee.

5. Conveying Modified Source Versions.

You may convey a work based on the Program, or the modifications to produce it from the Program, in the form of source code under the terms of section 4, provided that you also meet all of these conditions:

- (a) The work must carry prominent notices stating that you modified it, and giving a relevant date.
- (b) The work must carry prominent notices stating that it is released under this License and any conditions added under section 7. This requirement modifies the requirement in section 4 to "keep intact all notices".
- (c) You must license the entire work, as a whole, under this License to anyone who comes into possession of a copy. This License will therefore apply, along with any applicable section 7 additional terms, to the whole of the work, and all its parts, regardless of how they are packaged. This License gives no permission to license the work in any other way, but it does not invalidate such permission if you have separately received it.
- (d) If the work has interactive user interfaces, each must display Appropriate Legal Notices; however, if the Program has interactive interfaces that do not display Appropriate Legal Notices, your work need not make them do so.

A compilation of a covered work with other separate and independent works, which are not by their nature extensions of the covered work, and which are not combined with it such as to form a larger program, in or on a volume of a storage or distribution medium, is called an "aggregate" if the compilation and its resulting copyright are not used to limit the access or legal rights of the compilation's users beyond what the individual works permit. Inclusion of a covered work in an aggregate does not cause this License to apply to the other parts of the aggregate.

6. Conveying Non-Source Forms.

You may convey a covered work in object code form under the terms of sections 4 and 5, provided that you also convey the machine-readable Corresponding Source under the terms of this License, in one of these ways:

- (a) Convey the object code in, or embodied in, a physical product (including a physical distribution medium), accompanied by the Corresponding Source fixed on a durable physical medium customarily used for software interchange.

- (b) Convey the object code in, or embodied in, a physical product (including a physical distribution medium), accompanied by a written offer, valid for at least three years and valid for as long as you offer spare parts or customer support for that product model, to give anyone who possesses the object code either (1) a copy of the Corresponding Source for all the software in the product that is covered by this License, on a durable physical medium customarily used for software interchange, for a price no more than your reasonable cost of physically performing this conveying of source, or (2) access to copy the Corresponding Source from a network server at no charge.
- (c) Convey individual copies of the object code with a copy of the written offer to provide the Corresponding Source. This alternative is allowed only occasionally and noncommercially, and only if you received the object code with such an offer, in accord with subsection 6b.
- (d) Convey the object code by offering access from a designated place (gratis or for a charge), and offer equivalent access to the Corresponding Source in the same way through the same place at no further charge. You need not require recipients to copy the Corresponding Source along with the object code. If the place to copy the object code is a network server, the Corresponding Source may be on a different server (operated by you or a third party) that supports equivalent copying facilities, provided you maintain clear directions next to the object code saying where to find the Corresponding Source. Regardless of what server hosts the Corresponding Source, you remain obligated to ensure that it is available for as long as needed to satisfy these requirements.
- (e) Convey the object code using peer-to-peer transmission, provided you inform other peers where the object code and Corresponding Source of the work are being offered to the general public at no charge under subsection 6d.

A separable portion of the object code, whose source code is excluded from the Corresponding Source as a System Library, need not be included in conveying the object code work.

A “User Product” is either (1) a “consumer product”, which means any tangible personal property which is normally used for personal, family, or household purposes, or (2) anything designed or sold for incorporation into a dwelling. In determining whether a product is a consumer product, doubtful cases shall be resolved in favor of coverage. For a particular product received by a particular user, “normally used” refers to a typical or common use of that class of product, regardless of the status of the particular user or of the way in which the particular user actually uses, or expects or is expected to use, the product. A product is a consumer product regardless of whether the product has substantial commercial, industrial or non-consumer uses, unless such uses represent the only significant mode of use of the product.

“Installation Information” for a User Product means any methods, procedures, authorization keys, or other information required to install and execute modified versions of a covered work in that User Product from a modified version of its Corresponding Source. The information must suffice to ensure that the continued functioning of the modified object code is in no case prevented or interfered with solely because modification has been made.

If you convey an object code work under this section in, or with, or specifically for use in, a User Product, and the conveying occurs as part of a transaction in which the right of possession and use of the User Product is transferred to the recipient in perpetuity or for a fixed term (regardless of how the transaction is characterized), the Corresponding Source conveyed under this section must be accompanied by the Installation Information. But this requirement does not apply if neither you nor any third party retains the ability to install modified object code on the User Product (for example, the work has been installed in ROM).

The requirement to provide Installation Information does not include a requirement to continue to provide support service, warranty, or updates for a work that has been modified or installed by the recipient, or for the User Product in which it has been modified or installed. Access to a network may be denied when the modification itself materially and adversely affects the operation of the network or violates the rules and protocols for communication across the network.

Corresponding Source conveyed, and Installation Information provided, in accord with this section must be in a format that is publicly documented (and with an implementation available to the public in source code form), and must require no special password or key for unpacking, reading or copying.

7. Additional Terms.

“Additional permissions” are terms that supplement the terms of this License by making exceptions from one or more of its conditions. Additional permissions that are applicable to the entire Program shall be treated as though they were included in this License, to the extent that they are valid under applicable law. If additional permissions apply only to part of the Program, that part may be used separately under those permissions, but the entire Program remains governed by this License without regard to the additional permissions.

When you convey a copy of a covered work, you may at your option remove any additional permissions from that copy, or from any part of it. (Additional permissions may be written to require their own removal in certain cases when you modify the work.) You may place additional permissions on material, added by you to a covered work, for which you have or can give appropriate copyright permission.

Notwithstanding any other provision of this License, for material you add to a covered work, you may (if authorized by the copyright holders of that material) supplement the terms of this License with terms:

- (a) Disclaiming warranty or limiting liability differently from the terms of sections 15 and 16 of this License; or
- (b) Requiring preservation of specified reasonable legal notices or author attributions in that material or in the Appropriate Legal Notices displayed by works containing it; or
- (c) Prohibiting misrepresentation of the origin of that material, or requiring that modified versions of such material be marked in reasonable ways as different from the original version; or
- (d) Limiting the use for publicity purposes of names of licensors or authors of the material; or
- (e) Declining to grant rights under trademark law for use of some trade names, trademarks, or service marks; or
- (f) Requiring indemnification of licensors and authors of that material by anyone who conveys the material (or modified versions of it) with contractual assumptions of liability to the recipient, for any liability that these contractual assumptions directly impose on those licensors and authors.

All other non-permissive additional terms are considered “further restrictions” within the meaning of section 10. If the Program as you received it, or any part of it, contains a notice stating that it is governed by this License along with a term that is a further restriction, you may remove that term. If a license document contains a further restriction but permits relicensing or conveying under this License, you may add to a covered work material governed by the terms of that license document, provided that the further restriction does not survive such relicensing or conveying.

If you add terms to a covered work in accord with this section, you must place, in the relevant source files, a statement of the additional terms that apply to those files, or a notice indicating where to find the applicable terms.

Additional terms, permissive or non-permissive, may be stated in the form of a separately written license, or stated as exceptions; the above requirements apply either way.

8. Termination.

You may not propagate or modify a covered work except as expressly provided under this License. Any attempt otherwise to propagate or modify it is void, and will automatically terminate your rights under this License (including any patent licenses granted under the third paragraph of section 11).

However, if you cease all violation of this License, then your license from a particular copyright holder is reinstated (a) provisionally, unless and until the copyright holder explicitly and finally terminates your license, and (b) permanently, if the copyright holder fails to notify you of the violation by some reasonable means prior to 60 days after the cessation.

Moreover, your license from a particular copyright holder is reinstated permanently if the copyright holder notifies you of the violation by some reasonable means, this is the first time you have received

notice of violation of this License (for any work) from that copyright holder, and you cure the violation prior to 30 days after your receipt of the notice.

Termination of your rights under this section does not terminate the licenses of parties who have received copies or rights from you under this License. If your rights have been terminated and not permanently reinstated, you do not qualify to receive new licenses for the same material under section 10.

9. Acceptance Not Required for Having Copies.

You are not required to accept this License in order to receive or run a copy of the Program. Ancillary propagation of a covered work occurring solely as a consequence of using peer-to-peer transmission to receive a copy likewise does not require acceptance. However, nothing other than this License grants you permission to propagate or modify any covered work. These actions infringe copyright if you do not accept this License. Therefore, by modifying or propagating a covered work, you indicate your acceptance of this License to do so.

10. Automatic Licensing of Downstream Recipients.

Each time you convey a covered work, the recipient automatically receives a license from the original licensors, to run, modify and propagate that work, subject to this License. You are not responsible for enforcing compliance by third parties with this License.

An “entity transaction” is a transaction transferring control of an organization, or substantially all assets of one, or subdividing an organization, or merging organizations. If propagation of a covered work results from an entity transaction, each party to that transaction who receives a copy of the work also receives whatever licenses to the work the party’s predecessor in interest had or could give under the previous paragraph, plus a right to possession of the Corresponding Source of the work from the predecessor in interest, if the predecessor has it or can get it with reasonable efforts.

You may not impose any further restrictions on the exercise of the rights granted or affirmed under this License. For example, you may not impose a license fee, royalty, or other charge for exercise of rights granted under this License, and you may not initiate litigation (including a cross-claim or counterclaim in a lawsuit) alleging that any patent claim is infringed by making, using, selling, offering for sale, or importing the Program or any portion of it.

11. Patents.

A “contributor” is a copyright holder who authorizes use under this License of the Program or a work on which the Program is based. The work thus licensed is called the contributor’s “contributor version”.

A contributor’s “essential patent claims” are all patent claims owned or controlled by the contributor, whether already acquired or hereafter acquired, that would be infringed by some manner, permitted by this License, of making, using, or selling its contributor version, but do not include claims that would be infringed only as a consequence of further modification of the contributor version. For purposes of this definition, “control” includes the right to grant patent sublicenses in a manner consistent with the requirements of this License.

Each contributor grants you a non-exclusive, worldwide, royalty-free patent license under the contributor’s essential patent claims, to make, use, sell, offer for sale, import and otherwise run, modify and propagate the contents of its contributor version.

In the following three paragraphs, a “patent license” is any express agreement or commitment, however denominated, not to enforce a patent (such as an express permission to practice a patent or covenant not to sue for patent infringement). To “grant” such a patent license to a party means to make such an agreement or commitment not to enforce a patent against the party.

If you convey a covered work, knowingly relying on a patent license, and the Corresponding Source of the work is not available for anyone to copy, free of charge and under the terms of this License, through a publicly available network server or other readily accessible means, then you must either (1) cause the Corresponding Source to be so available, or (2) arrange to deprive yourself of the benefit of the

patent license for this particular work, or (3) arrange, in a manner consistent with the requirements of this License, to extend the patent license to downstream recipients. “Knowingly relying” means you have actual knowledge that, but for the patent license, your conveying the covered work in a country, or your recipient’s use of the covered work in a country, would infringe one or more identifiable patents in that country that you have reason to believe are valid.

If, pursuant to or in connection with a single transaction or arrangement, you convey, or propagate by procuring conveyance of, a covered work, and grant a patent license to some of the parties receiving the covered work authorizing them to use, propagate, modify or convey a specific copy of the covered work, then the patent license you grant is automatically extended to all recipients of the covered work and works based on it.

A patent license is “discriminatory” if it does not include within the scope of its coverage, prohibits the exercise of, or is conditioned on the non-exercise of one or more of the rights that are specifically granted under this License. You may not convey a covered work if you are a party to an arrangement with a third party that is in the business of distributing software, under which you make payment to the third party based on the extent of your activity of conveying the work, and under which the third party grants, to any of the parties who would receive the covered work from you, a discriminatory patent license (a) in connection with copies of the covered work conveyed by you (or copies made from those copies), or (b) primarily for and in connection with specific products or compilations that contain the covered work, unless you entered into that arrangement, or that patent license was granted, prior to 28 March 2007.

Nothing in this License shall be construed as excluding or limiting any implied license or other defenses to infringement that may otherwise be available to you under applicable patent law.

12. No Surrender of Others’ Freedom.

If conditions are imposed on you (whether by court order, agreement or otherwise) that contradict the conditions of this License, they do not excuse you from the conditions of this License. If you cannot convey a covered work so as to satisfy simultaneously your obligations under this License and any other pertinent obligations, then as a consequence you may not convey it at all. For example, if you agree to terms that obligate you to collect a royalty for further conveying from those to whom you convey the Program, the only way you could satisfy both those terms and this License would be to refrain entirely from conveying the Program.

13. Use with the GNU Affero General Public License.

Notwithstanding any other provision of this License, you have permission to link or combine any covered work with a work licensed under version 3 of the GNU Affero General Public License into a single combined work, and to convey the resulting work. The terms of this License will continue to apply to the part which is the covered work, but the special requirements of the GNU Affero General Public License, section 13, concerning interaction through a network will apply to the combination as such.

14. Revised Versions of this License.

The Free Software Foundation may publish revised and/or new versions of the GNU General Public License from time to time. Such new versions will be similar in spirit to the present version, but may differ in detail to address new problems or concerns.

Each version is given a distinguishing version number. If the Program specifies that a certain numbered version of the GNU General Public License “or any later version” applies to it, you have the option of following the terms and conditions either of that numbered version or of any later version published by the Free Software Foundation. If the Program does not specify a version number of the GNU General Public License, you may choose any version ever published by the Free Software Foundation.

If the Program specifies that a proxy can decide which future versions of the GNU General Public License can be used, that proxy’s public statement of acceptance of a version permanently authorizes you to choose that version for the Program.

Later license versions may give you additional or different permissions. However, no additional obligations are imposed on any author or copyright holder as a result of your choosing to follow a later version.

15. Disclaimer of Warranty.

THERE IS NO WARRANTY FOR THE PROGRAM, TO THE EXTENT PERMITTED BY APPLICABLE LAW. EXCEPT WHEN OTHERWISE STATED IN WRITING THE COPYRIGHT HOLDERS AND/OR OTHER PARTIES PROVIDE THE PROGRAM “AS IS” WITHOUT WARRANTY OF ANY KIND, EITHER EXPRESSED OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. THE ENTIRE RISK AS TO THE QUALITY AND PERFORMANCE OF THE PROGRAM IS WITH YOU. SHOULD THE PROGRAM PROVE DEFECTIVE, YOU ASSUME THE COST OF ALL NECESSARY SERVICING, REPAIR OR CORRECTION.

16. Limitation of Liability.

IN NO EVENT UNLESS REQUIRED BY APPLICABLE LAW OR AGREED TO IN WRITING WILL ANY COPYRIGHT HOLDER, OR ANY OTHER PARTY WHO MODIFIES AND/OR CONVEYS THE PROGRAM AS PERMITTED ABOVE, BE LIABLE TO YOU FOR DAMAGES, INCLUDING ANY GENERAL, SPECIAL, INCIDENTAL OR CONSEQUENTIAL DAMAGES ARISING OUT OF THE USE OR INABILITY TO USE THE PROGRAM (INCLUDING BUT NOT LIMITED TO LOSS OF DATA OR DATA BEING RENDERED INACCURATE OR LOSSES SUSTAINED BY YOU OR THIRD PARTIES OR A FAILURE OF THE PROGRAM TO OPERATE WITH ANY OTHER PROGRAMS), EVEN IF SUCH HOLDER OR OTHER PARTY HAS BEEN ADVISED OF THE POSSIBILITY OF SUCH DAMAGES.

17. Interpretation of Sections 15 and 16.

If the disclaimer of warranty and limitation of liability provided above cannot be given local legal effect according to their terms, reviewing courts shall apply local law that most closely approximates an absolute waiver of all civil liability in connection with the Program, unless a warranty or assumption of liability accompanies a copy of the Program in return for a fee.

END OF TERMS AND CONDITIONS

How to Apply These Terms to Your New Programs

If you develop a new program, and you want it to be of the greatest possible use to the public, the best way to achieve this is to make it free software which everyone can redistribute and change under these terms.

To do so, attach the following notices to the program. It is safest to attach them to the start of each source file to most effectively state the exclusion of warranty; and each file should have at least the “copyright” line and a pointer to where the full notice is found.

<one line to give the program’s name and a brief idea of what it does.>

Copyright (C) <textyear> <name of author>

This program is free software: you can redistribute it and/or modify it under the terms of the GNU General Public License as published by the Free Software Foundation, either version 3 of the License, or (at your option) any later version.

This program is distributed in the hope that it will be useful, but WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the GNU General Public License for more details.

You should have received a copy of the GNU General Public License along with this program. If not, see <<http://www.gnu.org/licenses/>>.

Also add information on how to contact you by electronic and paper mail.

If the program does terminal interaction, make it output a short notice like this when it starts in an interactive mode:

```
<program> Copyright (C) <year> <name of author>
```

```
This program comes with ABSOLUTELY NO WARRANTY; for details type 'show w'.  
This is free software, and you are welcome to redistribute it  
under certain conditions; type 'show c' for details.
```

The hypothetical commands `show w` and `show c` should show the appropriate parts of the General Public License. Of course, your program's commands might be different; for a GUI interface, you would use an “about box”.

You should also get your employer (if you work as a programmer) or school, if any, to sign a “copyright disclaimer” for the program, if necessary. For more information on this, and how to apply and follow the GNU GPL, see <http://www.gnu.org/licenses/>.

The GNU General Public License does not permit incorporating your program into proprietary programs. If your program is a subroutine library, you may consider it more useful to permit linking proprietary applications with the library. If this is what you want to do, use the GNU Lesser General Public License instead of this License. But first, please read <http://www.gnu.org/philosophy/why-not-lgpl.html>.

8.3 GNU Lesser General Public License v3

GNU LESSER GENERAL PUBLIC LICENSE

Version 3, 29 June 2007

Copyright © 2007 Free Software Foundation, Inc. <http://fsf.org/>

Everyone is permitted to copy and distribute verbatim copies of this license document, but changing it is not allowed.

This version of the GNU Lesser General Public License incorporates the terms and conditions of version 3 of the GNU General Public License, supplemented by the additional permissions listed below.

0. Additional Definitions.

As used herein, “this License” refers to version 3 of the GNU Lesser General Public License, and the “GNU GPL” refers to version 3 of the GNU General Public License.

“The Library” refers to a covered work governed by this License, other than an Application or a Combined Work as defined below.

An “Application” is any work that makes use of an interface provided by the Library, but which is not otherwise based on the Library. Defining a subclass of a class defined by the Library is deemed a mode of using an interface provided by the Library.

A “Combined Work” is a work produced by combining or linking an Application with the Library. The particular version of the Library with which the Combined Work was made is also called the “Linked Version”.

The “Minimal Corresponding Source” for a Combined Work means the Corresponding Source for the Combined Work, excluding any source code for portions of the Combined Work that, considered in isolation, are based on the Application, and not on the Linked Version.

The “Corresponding Application Code” for a Combined Work means the object code and/or source code for the Application, including any data and utility programs needed for reproducing the Combined Work from the Application, but excluding the System Libraries of the Combined Work.

1. Exception to Section 3 of the GNU GPL.

You may convey a covered work under sections 3 and 4 of this License without being bound by section 3 of the GNU GPL.

2. Conveying Modified Versions.

If you modify a copy of the Library, and, in your modifications, a facility refers to a function or data to be supplied by an Application that uses the facility (other than as an argument passed when the facility is invoked), then you may convey a copy of the modified version:

- (a) under this License, provided that you make a good faith effort to ensure that, in the event an Application does not supply the function or data, the facility still operates, and performs whatever part of its purpose remains meaningful, or
- (b) under the GNU GPL, with none of the additional permissions of this License applicable to that copy.

3. Object Code Incorporating Material from Library Header Files.

The object code form of an Application may incorporate material from a header file that is part of the Library. You may convey such object code under terms of your choice, provided that, if the incorporated material is not limited to numerical parameters, data structure layouts and accessors, or small macros, inline functions and templates (ten or fewer lines in length), you do both of the following:

- (a) Give prominent notice with each copy of the object code that the Library is used in it and that the Library and its use are covered by this License.

- (b) Accompany the object code with a copy of the GNU GPL and this license document.

4. Combined Works.

You may convey a Combined Work under terms of your choice that, taken together, effectively do not restrict modification of the portions of the Library contained in the Combined Work and reverse engineering for debugging such modifications, if you also do each of the following:

- (a) Give prominent notice with each copy of the Combined Work that the Library is used in it and that the Library and its use are covered by this License.
- (b) Accompany the Combined Work with a copy of the GNU GPL and this license document.
- (c) For a Combined Work that displays copyright notices during execution, include the copyright notice for the Library among these notices, as well as a reference directing the user to the copies of the GNU GPL and this license document.
- (d) Do one of the following:
 - . Convey the Minimal Corresponding Source under the terms of this License, and the Corresponding Application Code in a form suitable for, and under terms that permit, the user to recombine or relink the Application with a modified version of the Linked Version to produce a modified Combined Work, in the manner specified by section 6 of the GNU GPL for conveying Corresponding Source.
 - i. Use a suitable shared library mechanism for linking with the Library. A suitable mechanism is one that (a) uses at run time a copy of the Library already present on the user's computer system, and (b) will operate properly with a modified version of the Library that is interface-compatible with the Linked Version.
- (e) Provide Installation Information, but only if you would otherwise be required to provide such information under section 6 of the GNU GPL, and only to the extent that such information is necessary to install and execute a modified version of the Combined Work produced by recombining or relinking the Application with a modified version of the Linked Version. (If you use option 4d0, the Installation Information must accompany the Minimal Corresponding Source and Corresponding Application Code. If you use option 4d1, you must provide the Installation Information in the manner specified by section 6 of the GNU GPL for conveying Corresponding Source.)

5. Combined Libraries.

You may place library facilities that are a work based on the Library side by side in a single library together with other library facilities that are not Applications and are not covered by this License, and convey such a combined library under terms of your choice, if you do both of the following:

- (a) Accompany the combined library with a copy of the same work based on the Library, uncombined with any other library facilities, conveyed under the terms of this License.
- (b) Give prominent notice with the combined library that part of it is a work based on the Library, and explaining where to find the accompanying uncombined form of the same work.

6. Revised Versions of the GNU Lesser General Public License.

The Free Software Foundation may publish revised and/or new versions of the GNU Lesser General Public License from time to time. Such new versions will be similar in spirit to the present version, but may differ in detail to address new problems or concerns.

Each version is given a distinguishing version number. If the Library as you received it specifies that a certain numbered version of the GNU Lesser General Public License “or any later version” applies to it, you have the option of following the terms and conditions either of that published version or of any later version published by the Free Software Foundation. If the Library as you received it does not specify a version number of the GNU Lesser General Public License, you may choose any version of the GNU Lesser General Public License ever published by the Free Software Foundation.

If the Library as you received it specifies that a proxy can decide whether future versions of the GNU Lesser General Public License shall apply, that proxy's public statement of acceptance of any version is permanent authorization for you to choose that version for the Library.

8.4 MIT License

Permission is hereby granted, free of charge, to any person obtaining a copy of this software and associated documentation files (the "Software"), to deal in the Software without restriction, including without limitation the rights to use, copy, modify, merge, publish, distribute, sublicense, and/or sell copies of the Software, and to permit persons to whom the Software is furnished to do so, subject to the following conditions:

The above copyright notice and this permission notice shall be included in all copies or substantial portions of the Software.

THE SOFTWARE IS PROVIDED "AS IS", WITHOUT WARRANTY OF ANY KIND, EXPRESS OR IMPLIED, INCLUDING BUT NOT LIMITED TO THE WARRANTIES OF MERCHANTABILITY, FITNESS FOR A PARTICULAR PURPOSE AND NONINFRINGEMENT. IN NO EVENT SHALL THE AUTHORS OR COPYRIGHT HOLDERS BE LIABLE FOR ANY CLAIM, DAMAGES OR OTHER LIABILITY, WHETHER IN AN ACTION OF CONTRACT, TORT OR OTHERWISE, ARISING FROM, OUT OF OR IN CONNECTION WITH THE SOFTWARE OR THE USE OR OTHER DEALINGS IN THE SOFTWARE.

8.5 Zope Public License v2.1

A copyright notice accompanies this license document that identifies the copyright holders. This license has been certified as open source. It has also been designated as GPL compatible by the Free Software Foundation (FSF). Redistribution and use in source and binary forms, with or without modification, are permitted provided that the following conditions are met:

1. Redistributions in source code must retain the accompanying copyright notice, this list of conditions, and the following disclaimer.
2. Redistributions in binary form must reproduce the accompanying copyright notice, this list of conditions, and the following disclaimer in the documentation and/or other materials provided with the distribution.
3. Names of the copyright holders must not be used to endorse or promote products derived from this software without prior written permission from the copyright holders.
4. The right to distribute this software or to use it for any purpose does not give You the right to use Servicemarks (sm) or Trademarks (tm) of the copyright holders.
5. If You make a patent infringement claim against any entity (including a cross-claim or counterclaim in a lawsuit) alleging that the software licensed hereunder constitutes direct or contributory patent infringement, then any patent and any right-to-use licenses granted to You under this License for that software shall terminate as of the date such claim is made.

If any files are modified by You then you must cause the modified files to carry prominent notices stating that you changed the files and the date of any change.

Disclaimer

THIS SOFTWARE IS PROVIDED BY THE COPYRIGHT HOLDERS AS IS AND ANY EXPRESSED OR IMPLIED WARRANTIES, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE ARE DISCLAIMED. IN NO EVENT SHALL THE COPYRIGHT HOLDERS BE LIABLE FOR ANY DIRECT, INDIRECT, INCIDENTAL, SPECIAL, EXEMPLARY, OR CONSEQUENTIAL DAMAGES (INCLUDING, BUT NOT LIMITED TO, PROCUREMENT OF SUBSTITUTE GOODS

OR SERVICES; LOSS OF USE, DATA, OR PROFITS; OR BUSINESS INTERRUPTION) HOWEVER CAUSED AND ON ANY THEORY OF LIABILITY, WHETHER IN CONTRACT, STRICT LIABILITY, OR TORT (INCLUDING NEGLIGENCE OR OTHERWISE) ARISING IN ANY WAY OUT OF THE USE OF THIS SOFTWARE, EVEN IF ADVISED OF THE POSSIBILITY OF SUCH DAMAGE.

8.6 Apache License v2.0

Apache License
Version 2.0, January 2004
<http://www.apache.org/licenses/>

TERMS AND CONDITIONS FOR USE, REPRODUCTION, AND DISTRIBUTION

1. Definitions.

"License" shall mean the terms and conditions for use, reproduction, and distribution as defined by Sections 1 through 9 of this document.

"Licensor" shall mean the copyright owner or entity authorized by the copyright owner that is granting the License.

"Legal Entity" shall mean the union of the acting entity and all other entities that control, are controlled by, or are under common control with that entity. For the purposes of this definition, "control" means (i) the power, direct or indirect, to cause the direction or management of such entity, whether by contract or otherwise, or (ii) ownership of fifty percent (50) outstanding shares, or (iii) beneficial ownership of such entity.

"You" (or "Your") shall mean an individual or Legal Entity exercising permissions granted by this License.

"Source" form shall mean the preferred form for making modifications, including but not limited to software source code, documentation source, and configuration files.

"Object" form shall mean any form resulting from mechanical transformation or translation of a Source form, including but not limited to compiled object code, generated documentation, and conversions to other media types.

"Work" shall mean the work of authorship, whether in Source or Object form, made available under the License, as indicated by a copyright notice that is included in or attached to the work (an example is provided in the Appendix below).

"Derivative Works" shall mean any work, whether in Source or Object form, that is based on (or derived from) the Work and for which the editorial revisions, annotations, elaborations, or other modifications represent, as a whole, an original work of authorship. For the purposes of this License, Derivative Works shall not include works that remain separable from, or merely link (or bind by name) to the interfaces of, the Work and Derivative Works thereof.

"Contribution" shall mean any work of authorship, including the original version of the Work and any modifications or additions to that Work or Derivative Works thereof, that is intentionally submitted to Licensor for inclusion in the Work by the copyright owner or by an individual or Legal Entity authorized to submit on behalf of the copyright owner. For the purposes of this definition, "submitted" means any form of electronic, verbal, or written communication sent to the Licensor or its representatives, including but not limited to communication on electronic mailing lists, source code control systems, and issue tracking systems that are managed by, or on behalf of, the Licensor for the purpose of discussing and improving the Work, but excluding communication that is conspicuously marked or otherwise designated in writing by the copyright owner as "Not a Contribution."

"Contributor" shall mean Licensor and any individual or Legal Entity on behalf of whom a Contribution has been received by Licensor and subsequently incorporated within the Work.

2. Grant of Copyright License. Subject to the terms and conditions of this License, each Contributor hereby grants to You a perpetual, worldwide, non-exclusive, no-charge, royalty-free, irrevocable copyright license to reproduce, prepare Derivative Works of, publicly display, publicly perform, sublicense, and distribute the Work and such Derivative Works in Source or Object form.

3. Grant of Patent License. Subject to the terms and conditions of this License, each Contributor hereby grants to You a perpetual, worldwide, non-exclusive, no-charge, royalty-free, irrevocable (except as stated in this section) patent license to make, have made, use, offer to sell, sell, import, and otherwise transfer the Work, where such license applies only to those patent claims licensable by such Contributor that are necessarily infringed by their Contribution(s) alone or by combination of their Contribution(s) with the Work to which such Contribution(s) was submitted. If You institute patent litigation against any entity (including a cross-claim or counterclaim in a lawsuit) alleging that the Work or a Contribution incorporated within the Work constitutes direct or contributory patent infringement, then any patent licenses granted to You under this License for that Work shall terminate as of the date such litigation is filed.

4. Redistribution. You may reproduce and distribute copies of the Work or Derivative Works thereof in any medium, with or without modifications, and in Source or Object form, provided that You meet the following conditions:

- (a) You must give any other recipients of the Work or Derivative Works a copy of this License; and
- (b) You must cause any modified files to carry prominent notices stating that You changed the files; and
- (c) You must retain, in the Source form of any Derivative Works that You distribute, all copyright, patent, trademark, and attribution notices from the Source form of the Work, excluding those notices that do not pertain to any part of the Derivative Works; and

- (d) If the Work includes a "NOTICE" text file as part of its distribution, then any Derivative Works that You distribute must include a readable copy of the attribution notices contained within such NOTICE file, excluding those notices that do not pertain to any part of the Derivative Works, in at least one of the following places: within a NOTICE text file distributed as part of the Derivative Works; within the Source form or documentation, if provided along with the Derivative Works; or, within a display generated by the Derivative Works, if and wherever such third-party notices normally appear. The contents of the NOTICE file are for informational purposes only and do not modify the License. You may add Your own attribution notices within Derivative Works that You distribute, alongside or as an addendum to the NOTICE text from the Work, provided that such additional attribution notices cannot be construed as modifying the License.

You may add Your own copyright statement to Your modifications and may provide additional or different license terms and conditions for use, reproduction, or distribution of Your modifications, or for any such Derivative Works as a whole, provided Your use, reproduction, and distribution of the Work otherwise complies with the conditions stated in this License.

5. Submission of Contributions. Unless You explicitly state otherwise, any Contribution intentionally submitted for inclusion in the Work by You to the Licensor shall be under the terms and conditions of this License, without any additional terms or conditions. Notwithstanding the above, nothing herein shall supersede or modify the terms of any separate license agreement you may have executed with Licensor regarding such Contributions.

6. Trademarks. This License does not grant permission to use the trade names, trademarks, service marks, or product names of the Licensor, except as required for reasonable and customary use in describing the origin of the Work and reproducing the content of the NOTICE file.

7. Disclaimer of Warranty. Unless required by applicable law or agreed to in writing, Licensor provides the Work (and each Contributor provides its Contributions) on an "AS IS" BASIS, WITHOUT WARRANTIES OR CONDITIONS OF ANY KIND, either express or implied, including, without limitation, any warranties or conditions of TITLE, NON-INFRINGEMENT, MERCHANTABILITY, or FITNESS FOR A PARTICULAR PURPOSE. You are solely responsible for determining the appropriateness of using or redistributing the Work and assume any risks associated with Your exercise of permissions under this License.

8. Limitation of Liability. In no event and under no legal theory, whether in tort (including negligence), contract, or otherwise, unless required by applicable law (such as deliberate and grossly negligent acts) or agreed to in writing, shall any Contributor be liable to You for damages, including any direct, indirect, special, incidental, or consequential damages of any character arising as a result of this License or out of the use or inability to use the Work (including but not limited to damages for loss of goodwill, work stoppage, computer failure or malfunction, or any and all other commercial damages or losses), even if such Contributor has been advised of the possibility of such damages.

9. Accepting Warranty or Additional Liability. While redistributing the Work or Derivative Works thereof, You may choose to offer, and charge a fee for, acceptance of support, warranty, indemnity, or other liability obligations and/or rights consistent with this License. However, in accepting such obligations, You may act only on Your own behalf and on Your sole responsibility, not on behalf of any other Contributor, and only if You agree to indemnify, defend, and hold each Contributor harmless for any liability incurred by, or claims asserted against, such Contributor by reason of your accepting any such warranty or additional liability.

END OF TERMS AND CONDITIONS

APPENDIX: How to apply the Apache License to your work.

To apply the Apache License to your work, attach the following boilerplate notice, with the fields enclosed by brackets "[]" replaced with your own identifying information. (Don't include the brackets!) The text should be enclosed in the appropriate comment syntax for the file format. We also recommend that a file or

class name and description of purpose be included on the same "printed page" as the copyright notice for easier identification within third-party archives.

Copyright [yyyy] [name of copyright owner]

Licensed under the Apache License, Version 2.0 (the "License");
you may not use this file except in compliance with the License.
You may obtain a copy of the License at

<http://www.apache.org/licenses/LICENSE-2.0>

Unless required by applicable law or agreed to in writing, software
distributed under the License is distributed on an "AS IS" BASIS,
WITHOUT WARRANTIES OR CONDITIONS OF ANY KIND, either express or implied.
See the License for the specific language governing permissions and
limitations under the License.